EOS-X

BY DĆZSSZ PRO™

Operating Manual

It is mandatory that users or operators of this device read and comprehend all information in this operating manual prior to using the EOS-X by Déesse PRO device.

You must follow all the instructions in this operating manual.

Please contact the brand owners, Harpar Grace International, if you have any comments about the EOS-X by Déesse PRO device or this operating manual.

Please contact Harpar Grace International if patients experience any side effects that are not mentioned or are more severe than those described in the operating manual.

☎ +44 2038 686242 ✉ enquiries@harpargrace.com

Harpar Grace International Limited is the brand owner of Déesse PRO. Déesse PRO is the trademark of Harpar Grace International Limited.

Harpar Grace International Ltd, Unit 4 Imperial Court, Magellan Close Walworth Business Park Andover, Hampshire, SP10 5NT, United Kingdom

A copy of this operating manual is available at www.deessepro.com

Certainty IP Ltd, Dublin, D18 CV48. +353 1582 9003

Emergo Australia, Level 20 Tower II, Darling Park, 201 Sussex Street, Sydney, NSW 2000, Australia

Legal Manufacturer:

iSmart Developments Ltd, 129 Green Lanes, Sutton Coldfield, West Midlands, B73 5LT, United Kingdom

The information in this document is proprietary to iSmart Developments Ltd, and supplied on the understanding that it shall be kept confidential to recipients and that it shall not be reproduced and/or used for commercial purposes without the written consent of iSmart Developments Ltd.
21st November 2023 DOC 4281 V0.

Contents

Contents	3
1.0. Introduction	4
2.0. Device Description	4
3.0. Intended Use	5
4.0. Contraindications, Warnings & Precautions	6
4.1. Contraindications	6
4.2. Warnings	7
4.3. General Safety	8
4.4. Precautions	8
5.0. Known Side Effects	9
6.0. Clinical Benefits	9
7.0. First-Time Setup	10
7.1. Unpacking the EOS-X by Déesse PRO	10
7.2. Connecting the Déesse PRO Face Mask & Neck Attachments to the Controller	10
7.3. Understanding the Controller Screen	10
8.0. Using the EOS-X by Déesse PRO Device	12
8.1. Treatment Recommendations	12
8.2. Positioning the EOS-X by Déesse PRO Device	12
8.3. Carrying Out a Treatment	12
8.4. Treatment Modes & Wavelengths	13
9.0. Care & Maintenance	14
10.0. Cleaning	14
11.0. Disposal & Recycling	15
12.0. Troubleshooting	15
13.0. Labelling	16
14.0. Technical Specification	16
15.0. Explanation of Symbols	17
16.0. Warranty	18

1.0. Introduction


Déesse PRO is the frontier in advanced light therapy, pioneering the intersection between clinical and at-home experiences. We combine over twenty years of clinical research and NASA-accredited principles to deliver purposeful innovation. Developed with a profound respect for balance, simplicity and renewal, our treatments optimise the skin and restore the mind to an essential state of equilibrium.

Clinically proven to stimulate cellular regeneration, our FDA-cleared light technology targets multiple skin and wellbeing concerns, resulting in both immediate and long-term benefits to support a lifetime of needs. Bespoke treatment programs and continual expert guidance empower users to experience the spectrum of Déesse PRO's capabilities; delivering transformative skin results and enhanced user wellbeing in a results-driven and purposeful ritual.

Thank you for choosing Déesse PRO.

2.0. Device Description

1. Hard shell face mask with integrated eye protection
2. Hard shell neck attachment
3. Touchscreen controller unit
4. Adjustable mask straps
5. Adjustable neck straps
6. Mains power supply with country-specific adapters
7. Display stand
8. Operating manual
9. Micro SD Card

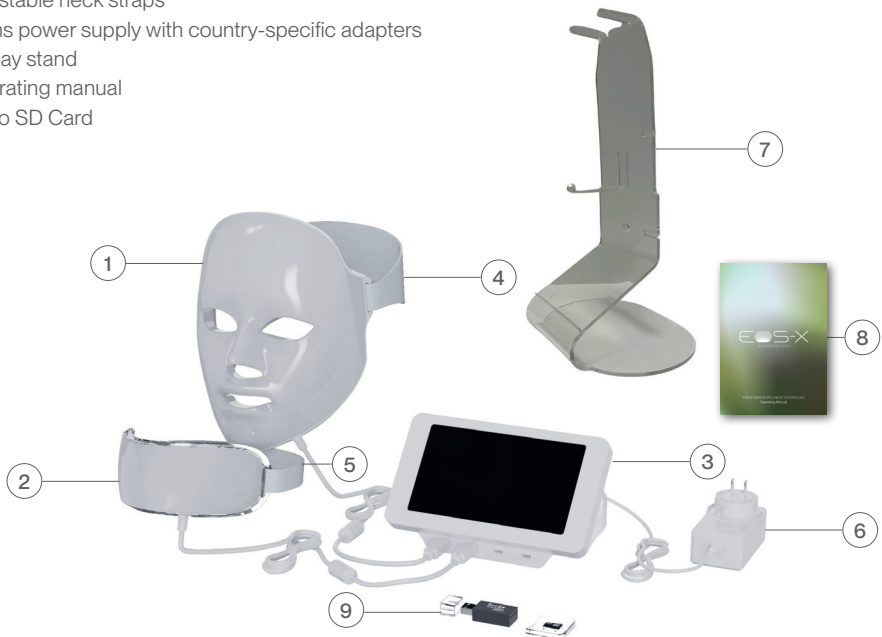Déesse PRO is a professional-grade LED light therapy device for treatment of the face and neck.

This professional device is intended for use in clinic, spa and salon environments and has multiple functions and protocols.

The device consists of a face mask containing light-emitting diodes (LEDs) ①, a neck attachment containing light-emitting diodes (LEDs) ② and a touchscreen controller unit ③. The device is powered by a 100-240V power adaptor ⑥.

The face mask is worn on the face and is held in place by adjustable straps ④. The neck attachment is worn on the neck and is held in place by the adjustable straps ⑤.

The touchscreen controller unit supplies power to the face mask and neck attachment and switches the LEDs ON/OFF and allows the operator to select the appropriate treatment modes. The face mask and neck attachment connect to the controller unit by separate USB C connectors. The power cable connects to the rear of the controller unit and the plug is plugged into a mains power supply. The country-specific adapters connect to the main plug for the relevant country-specific power outlet.

The face mask has 7 modes and emits light in the blue, red, green, and near-infrared regions of the light spectrum.

The neck attachment has 2 modes and emits light in the blue, red, and near-infrared regions of the light spectrum.

Should there be a requirement for a software update in the future, the device is supplied with a blank SD card. The SD card reader can also be used to update the controller with the release of any new LED device attachments in the future.

Please note, the equipment is not used to make measurements of any sort or to draw any conclusions regarding the indication to treat. The equipment does not require checks on the light output as the LEDs do not dim with age to any practical extent.

3.0. Intended Use

The EOS-X by Déesse PRO is a professional use light emitting diode device that is intended to emit light in the blue (415nm) and red (630nm) regions of the light spectrum for the treatment of mild to moderate acne vulgaris.

The EOS-X by Déesse PRO device is intended to emit light in the red (630nm) and near infra-red (NIR) region (830nm) of the light spectrum for the treatment of facial wrinkles.

Including the wavelengths above, the EOS-X by Déesse PRO device is also intended to emit light in the green (545nm) and near infra-red (1070nm) region of the light spectrum for the treatment of several

cosmetic indications and wellness benefits. This includes helping to improve the visible signs of skin ageing, improve skin tone and texture, reduce brown spots, and reduce inflammation and skin redness. The EOS-X by Déesse PRO device is also intended to help boost serotonin release, improve sleep quality, and help improve circadian rhythm.

4.0. Contraindications, Warnings & Precautions

4.1. Contraindications

- DO NOT use the EOS-X by Déesse PRO device to treat any other medical conditions apart from those listed in the intended use (full facial wrinkles and mild to moderate acne vulgaris). The EOS-X by Déesse PRO device has not been tested for any other conditions than those listed and the risk is unknown.
- DO NOT use the EOS-X by Déesse PRO device on women who are pregnant, trying to conceive, or breastfeeding. The EOS-X by Déesse PRO device has NOT been tested on pregnant or breast-feeding women and therefore the risk to the foetus, neonate or pregnant woman is unknown.
- DO NOT use the EOS-X by Déesse PRO device on anyone who has an allergy to polyester, nylon and Lycra. If they are allergic to any of these materials, you may not be able to use the mask and neck straps provided.
- DO NOT use the EOS-X by Déesse PRO device face mask and neck attachment on anyone who is allergic to Polycarbonate and/or ABS Plastics. If they have an allergy to these materials, you will not be able to use the EOS-X by Déesse PRO device.
- DO NOT use the EOS-X by Déesse PRO device on children under 14 years of age. The safety of the device for children younger than 14 years of age has NOT been tested and the risk is unknown.
- DO NOT use the EOS-X by Déesse PRO device on open cuts, wounds, oozing acne pustules and broken skin including herpes simplex infections.
- DO NOT use the EOS-X by Déesse PRO device on anyone who suffers from light-induced headaches.
- DO NOT use the EOS-X by Déesse PRO device on anyone who suffers from genetic conditions of the eyes. The EOS-X by Déesse PRO device has NOT been tested on subjects who suffer from genetic conditions of the eyes and therefore the risks are unknown. If you are unsure about any related eye condition and using the EOS-X by Déesse PRO device, contact a healthcare professional.
- DO NOT use the EOS-X by Déesse PRO device on anyone who suffers from any photosensitive disorder (sensitisation to light) of the eye. The EOS-X by Déesse PRO device has NOT been tested on subjects who suffer with photosensitive disorders of the eye and therefore the risks are unknown. If you are unsure about any related eye condition and using the EOS-X by Déesse PRO device, contact a healthcare professional.
- DO NOT use the EOS-X by Déesse PRO device on anyone who suffers from lupus erythematosus, photosensitive eczema or albinism. If you use the EOS-X by Déesse PRO device to treat someone who suffers from lupus erythematosus, photosensitive eczema or albinism, this may cause a severe skin reaction.
- DO NOT use the EOS-X by Déesse PRO device on anyone who suffers from any photosensitive disorder (sensitisation to light) of the skin. If you use the EOS-X by Déesse PRO device to treat someone who suffers from a photosensitive disorder, this may cause a severe reaction.

- DO NOT use the EOS-X by Déesse PRO device on anyone who is taking any medication that can cause photosensitivity. If you use the EOS-X by Déesse PRO device to treat someone who is taking any medication that can cause photosensitivity, this may cause a severe skin reaction.
- All practitioners must complete a treatment consent form with their clients before administering a treatment using the EOS-X by Déesse PRO device. A consent form template can be found in section 8.3.

For a full list of drugs and other substances that can cause photosensitivity, please visit: www.deessepro.com.

Photosensitivity is a common side effect of various medications. These can include certain antibiotics, chemotherapy drugs, and diuretics. If you are unsure about any medication you may be taking, consult your healthcare provider.

Other substances not listed above can also cause photosensitivity.

Common examples of these substances are St John's Wort, coal tar, deodorants, antibacterial soaps, artificial sweeteners, naphthalene (mothballs), petroleum products, brightening agents found in laundry detergent, and cadmium sulphide (a chemical injected into the skin during tattooing).

4.2. Warnings

- DO NOT use the EOS-X by Déesse PRO device in any way that is not exactly as specified in this operating manual. Use of the EOS-X by Déesse PRO device other than exactly as specified in this operating manual has not been tested and the risk is unknown.
- DO NOT use the EOS-X by Déesse PRO device on acne pustules that are exuding discharge or bleeding. Stop treatment and resume when the pustules have formed a plug or scab.
- DO NOT use the EOS-X by Déesse PRO device in an environment that is NOT specified by the manufacturer, as this may cause damage to the device or harm to the operator or patient.
- TURN OFF and REMOVE the EOS-X by Déesse PRO device immediately if the device or the skin feels hot. DO NOT continue to use the device and return to your place of purchase.
- STOP using the EOS-X by Déesse PRO device during a treatment if the patient complains of severe headache, changes in vision or pain in their eyes. Should this occur, contact a healthcare professional.
- DO NOT use the EOS-X by Déesse PRO device without first removing sunscreens, topical creams, lotions, or make-up BEFORE treatment. Some ingredients in these products may cause light sensitisation or may reduce the effect of the treatment. If you use the EOS-X by Déesse PRO device without removing these products, this may cause a skin reaction or may reduce the effect of the treatment.
- NEVER use the EOS-X by Déesse PRO device near water or where water can be spilt onto the equipment. DO NOT use the EOS-X by Déesse PRO device if you think it has come into contact with water.
- Keep the EOS-X by Déesse PRO device away from hot surfaces and DO NOT store or use the device in direct sunlight.
- DO NOT use the EOS-X by Déesse PRO device if it has been exposed to environmental conditions outside of those specified in the Technical Specification (see section 14.0.).

- DO NOT attempt to modify the EOS-X by Déesse PRO device in any way. Modifying the EOS-X by Déesse PRO device will void your warranty and can cause an electric shock when the system is used again. The EOS-X by Déesse PRO device system contains no serviceable parts.
- DO NOT use the EOS-X by Déesse PRO device if the face mask, neck attachment, controller unit, power plug, cables or packaging are damaged in any way.
- DO NOT use the EOS-X by Déesse PRO device if it shows any signs of damage to the material through wear and tear. If you follow the care and maintenance for your device, then the EOS-X by Déesse PRO device has a minimum lifetime of 2 years.
- DO NOT dispose of any part of the EOS-X by Déesse PRO device with normal household waste as this could cause environmental damage from the components contained in the device.

4.3. General Safety

The EOS-X by Déesse PRO device is an electrical device, which under certain circumstances could present an electrical shock hazard. Please read this operating manual thoroughly and follow the directions stated in the manual to assure maximum safety during operation.

EOS-X by Déesse PRO device may only be used with the power supply that is supplied with the device. The use of other power supplies may not be safe and could lead to electric shock.

DO NOT place the device in a location which makes it difficult to operate the device.

If you need to urgently terminate the treatment, you should use the on/off switch on the rear of the controller unit. This will instantly terminate the treatment.

4.4. Precautions

- DO NOT stare directly at the LED light source.
- DO NOT pick or rub acne pustules that have scabbed or blocked. Picking acne pustules may lead to acne scarring or worsening.
- Before and after each use, with the device switched off, clean the face mask and neck attachment thoroughly and properly. All surfaces of the EOS-X by Déesse PRO device must be cleaned by following the cleaning instructions in section 10.0., to prevent cross-contamination.
- DO NOT use any other cleaning products other than those stated in the cleaning instructions. Use of other cleansing products may damage the EOS-X by Déesse PRO device.
- Store the EOS-X by Déesse PRO device at room temperature in dry conditions, away from direct sunlight.
- The surface temperature of the device may exceed 41°C under maximum environmental operating temperature conditions
- Take care when using the EOS-X by Déesse PRO device around young children. Electrical cables can present a strangulation hazard and there is a risk of fire and burns due to children mishandling the EOS-X by Déesse PRO device.

- Check if the voltage indicated on the power supply provided corresponds to the local mains voltage before you connect the EOS-X by Déesse PRO device to avoid the risk of electrocution or permanent damage to the device.
- Make sure you plug the EOS-X by Déesse PRO device into a properly installed wall socket. If the wall socket used to power the EOS-X by Déesse PRO device has poor connections or the plug of the equipment becomes hot, there is a risk of fire and burns.
- DO NOT leave the EOS-X by Déesse PRO device unattended when it is switched on to avoid the risk of fire or burns.
- DO NOT subject the EOS-X by Déesse PRO device to heavy shocks or vibrations to avoid the risk of damage to the device. If you drop and damage the device, DO NOT use it any longer.
- DO NOT use the EOS-X by Déesse PRO device near other electromagnetic compatibility (EMC), electromagnetic fields (EMF) or radio frequency (RF) equipment, such as pacemakers to avoid potential interference and disturbances. These interferences and disturbances could cause malfunction of essential medical equipment.
- DO NOT use the device to overtreat - only use per the recommend treatment protocols

5.0. Known Side Effects

- Immediately after treatment skin may appear red. The skin should return to normal within 1-2 hours.
- Some users have reported pigmentation, skin flaking or peeling, associated pruritis and dryness after treatment of acne vulgaris. This is transient and clears within a few days.
- Temporary headaches immediately following the light treatment.

If your patient experiences a side effect that is greater than that described above or a serious side effect not listed, stop using the device immediately and seek advice from a healthcare provider. Report such an incident firstly to your country distributor (a list of authorised Déesse PRO distributors can be found on www.deessepro.com) or your purchase provider and your representative Competent Authority if you deem it applicable.

Although this device emits levels of visible red, blue, green and near-infrared light that are considered to be safe for the skin, if an adverse incident occurs after using the device, you should stop treatment immediately.

6.0. Clinical Benefits

People who have undergone a course of 2 weekly treatments for 5 weeks using the EOS-X by Déesse PRO device report a reduction of mild to moderate acne vulgaris.

People who have undergone a course of 15 treatments over a period of 6 weeks using the EOS-X by Déesse PRO device report a reduction in wrinkles.

7.0. First-Time Setup

7.1. Unpacking the EOS-X by Déesse PRO

Unpack the EOS-X by Déesse PRO device from the box and check off the parts described in the device description (section 2.0.).

Visually inspect the parts and connections. Return the device to your supplier if you think the device or connecting parts or components are damaged or have been previously used.

DO NOT use the EOS-X by Déesse PRO device if you notice the device, packaging, connecting parts or components are damaged.

7.2. Connecting the Déesse PRO Face Mask & Neck Attachments to the Controller

- Select the correct country adaptor for your country and connect by inserting the adaptor into the rear of the adaptor body and rotating clockwise until you hear an audible click. To remove, press down on the button and rotate anti-clockwise.
- Connect the power adaptor to the rear of the control panel via the coaxial jack connector.
- Connect the plug end of the power adaptor to a suitable wall socket.
- Switch ON the control panel by pressing the power button on the rear of the device.
- The control panel will display the EOS-X by Déesse PRO screen. Tap the screen to proceed to the main home screen.
- Connect the EOS-X by Déesse PRO face mask by inserting the USB C connector of the face mask into the USB C port on the front of the controller.
- Connect the EOS-X by Déesse PRO neck attachment by inserting the USB C connector of the neck attachment into the USB C port on the front of the controller.
- Select the cog settings icon in the menu to be able to choose your preferred language, volume and screen brightness.

7.3. Understanding the Controller Screen

Turn on the EOS-X by Déesse PRO device following section 7.2.

The touch screen contains options to allow navigation, change of settings, treatment selection and power control.

Menu Buttons

Use the following as a key to the function of the touchscreen buttons in the menu bar.

Home – this button returns the user to the mode selection screen

Information – this button gives the user access to warranty registration and training applications

Back – this button returns the user to the previous screen

Settings – this button takes the user to the settings screen where language, volume and screen brightness can be set and changed.

Sleep – this button will give the user the option to power down the controller into sleep mode

Treatment Screen Buttons

Use the following as a key to the function of the touchscreen buttons in the treatment screens.

Information for Treatment Additions – this button gives the user further information about the additional 1070nm booster treatment and the neck treatment options

Toggle – this button allows the user to enable or disable the additional treatment options

Information for User – this button gives the user access to the latest contraindications, warnings, precautions and other information via a QR code

Start – this button allows the user to start the treatment

Cancel – this button allows the user to cancel the treatment in progress

Pause – this button allows the user to pause the treatment in progress

Changing the Settings

The language, volume and screen brightness of the EOS-X by Déesse PRO controller can be changed in the settings screen by pressing the settings icon on the menu bar.

- **Language:** The touchscreen controller can be set in different languages. Choose your language by tapping on the desired option.
- **Volume:** The controller has an optional sound output to enhance the user experience, such as when the controller is powered on and when a treatment is complete. The volume of this sound output can be adjusted. Tap and drag the slider bar right to increase the volume and drag the slide left to decrease or mute the volume.
- **Brightness:** The brightness of the touchscreen can be adjusted by dragging the slider bar to the right to increase the brightness and dragging the slider left to dim the brightness.

8.0. Using the EOS-X by Déesse PRO Device

Regarding electromagnetic compatibility (EMC) and electromagnetic fields (EMF), in and around the use environment the Déesse PRO should be kept away from radio frequency (RF) equipment to avoid interference.

8.1. Treatment recommendations

For the treatment of mild to moderate acne vulgaris the best results are seen by treating for 2 times a week for 5 weeks.

For the treatment of wrinkles, the best results are seen by treating for 15 sessions over a period of 6 weeks.

Maintenance treatment should be carried out as required, according to lifestyle and skin condition.

8.2. Positioning the EOS-X by Déesse PRO Device

- The EOS-X by Déesse PRO device has two-part adjustable straps that allow you to position the face mask and neck attachment comfortably.
- To fit the face mask, thread the straps through the inlet holes on either side of the face mask. Loop the straps around the head and secure firmly with the Velcro at the rear.
- To fit the neck attachment, thread the straps through the inlet holes on either side of the neck attachment. Loop the straps around the neck and secure firmly with the Velcro at the rear.
- Adjust the straps if necessary for comfort.

8.3. Carrying Out a Treatment

All practitioners must complete a treatment consent form with their clients before administering a treatment using the EOS-X by Déesse PRO device. A consent form template can be downloaded by scanning the QR code.

For returning patients, ensure that you have reviewed any changes in their health or medication in line with their consent form.

The client is not intended to be the operator of this device.

1. Cleanse and exfoliate the skin to remove any traces of sunscreen, skincare, lotions or topical creams.
2. Fit the EOS-X by Déesse PRO face mask and/or neck attachment on the face or neck, and secure using the adjustable straps.
3. Connect the face mask and/or neck attachment to the controller unit and turn on the power button on the rear of the controller.
4. From the home screen on the controller unit, follow the instructions to select the treatment mode.

5. Selecting the mode will take you to the pre-treatment screen which provides you with the option of selecting treatment additions using the toggle buttons.
6. Once selected, press the start icon and the treatment will commence after a 3-second countdown.
7. The treatment screen displays the treatment mode and the time remaining.
8. The treatment can be paused or stopped at any time using the pause or cancel buttons on the screen. You can also end the treatment at any time by tapping the sleep icon.
9. At the end of the treatment the LEDs will switch OFF and a treatment complete screen will be displayed with an option to return to the home screen where sleep mode can be activated.
10. The controller screen will revert to the home screen automatically after 5 minutes.
11. Remove the face mask and/or neck attachment. The user may apply normal skincare or topicals as usual.
12. Unplug and clean the device, following the cleaning instructions in section 10.0.

8.4. Treatment Modes & Wavelengths

All treatment modes are set at the optimum treatment time of 10 minutes.

For further treatment guidance and protocols, please visit www.deessepro.com or scan the QR code.

Mask Treatment Modes	Medical Indication	Wavelengths	Treatment Additions
Mode 1: Anti-Ageing	Treatment of Facial Wrinkles (See Intended Use Statement)* (* USA only)	630nm (Red) 830nm (Near-Infrared)	Boost with 1070nm (Near – Infrared) Add Anti-Ageing Neck Treatment
Mode 2: Acne	Treatment of mild to moderate acne vulgaris. (See Intended Use Statement)	415nm (Blue) 630nm (Red)	Add Spots & Blemishes Neck Treatment
Mode 3: Brightening	N/A	545nm (Green)	Boost with 1070nm (Near – Infrared) Add Anti-Ageing Neck Treatment
Mode 4: Photoaging	N/A	630nm (Red) 830nm (Near-Infrared)	Boost with 1070nm (Near – Infrared) Add Anti-Ageing Neck Treatment
Mode 5: Skin Calming	N/A	415nm (Blue) 630nm (Red) 830nm (Near-Infrared)	Add Anti-Ageing Neck Treatment
Mode 6: Mood Boosting	N/A	415nm (Blue) 630nm (Red)	Add Anti-Ageing Neck Treatment
Mode 7: Relaxing	N/A	630nm (Red)	Add Anti-Ageing Neck Treatment

Modes 1, 3 and 4 are suitable to be boosted with the addition of 1070nm wavelength. This wavelength helps:

- Stimulate deeper layers of the skin
- Reduces fine line & wrinkles around eyes
- Reduce eye-bags & dark circles
- Boost the immune system

Neck Treatment Modes	Wavelengths
Mode 8: Anti-Ageing	630nm (Red) 830nm (Near-Infrared)
Mode 9: Blemishes & Spots	415nm (Blue) 630nm (Red)

9.0. Care & Maintenance

- Store the EOS-X by Déesse PRO device flat and at room temperature in dry conditions. Use the display stand or gift box provided to store your device when not in use.
- Store and use the EOS-X by Déesse PRO device in conditions set out in the labelling provided (see section 13.0.). Avoid storing and using the EOS-X by Déesse PRO device in dusty conditions, in direct sunlight or in damp and humid conditions. DO NOT use the EOS-X by Déesse PRO device if it has been exposed to environmental conditions outside of those specified in the Technical Specification (see section 14.0.).
- The EOS-X by Déesse PRO device system contains no serviceable parts. DO NOT attempt to open and try and modify the EOS-X by Déesse PRO device in any way. Hazards, such as electrical shock or burns, can result from unauthorised modification of the device. Modifying the EOS-X by Déesse PRO device will void your warranty.
- Before each use, with the device turned off, clean the device following the cleaning instructions (see section 10.0.).
- Before each use, visually inspect the EOS-X by Déesse PRO device face mask, neck attachment, controller unit and cables for any damage. DO NOT use the device if you notice any signs of visible damage.
- DO NOT use the EOS-X by Déesse PRO device if it shows any signs of damage to the material through wear and tear. If you follow the care and maintenance for your EOS-X by Déesse PRO, then the device has a minimum lifetime of 2 years.
- DO NOT subject the EOS-X by Déesse PRO device to heavy shocks or vibrations to avoid the risk of damage to the device. If you drop and damage the device, DO NOT use it any longer.

10.0. Cleaning

DO NOT use any other cleaning products on your EOS-X by Déesse PRO device as these may damage the materials.

Before cleaning ensure that the EOS-X by Déesse PRO device is unplugged from the mains power outlet.

Before and after use of the EOS-X by Déesse PRO, clean all surfaces of the face mask, neck attachment and controller unit thoroughly and properly with a sanitising/disinfecting wipe.

Following the cleaning procedure in this operating manual will NOT cause loss of basic safety performance.

11.0. Disposal & Recycling

The EOS-X by Déesse PRO device contains electrical components that are not suitable for disposal in household waste. DO NOT dispose of the EOS-X by Déesse PRO device in household waste.

Follow your country's rules for the separate collection of electrical and electronic products. Correct disposal helps prevent damage to the environment and human health. Please contact your country distributor (a list of authorised EOS-X by Déesse PRO device distributors can be found on www.deessepro.com) if you would like further information on the proper disposal of your EOS-X by Déesse PRO device.

The outer box packaging and this operating manual are recyclable. The device components (face mask, neck attachment, straps, controller unit, power supply and connecting cables) are not and should not be disposed of in household waste

12.0. Troubleshooting

Problem	Possible Cause	Action
Device does not turn on	Faulty power adaptor	Remove power adaptor from power outlet and visually inspect the power adaptor and cables for damage.
Device does not turn on	Faulty electrical input	If the power adaptor is NOT damaged ensure that the power adaptor is inserted into a functioning power socket.
Device does not turn on	Poor connection	If the power adaptor is NOT damaged ensure that the coaxial jack connector is pushed firmly into the control panel
User interface does not respond or becomes corrupted	Touchscreen is faulty	Power down the controller and unplug from the power outlet. Wait 30 seconds. Reconnect and turn on the controller.

Problem	Possible Cause	Action
The controller fails to register or recognise an attachment	Loose connection/faulty cable	Try removing and re-attaching the device cable into the controller USB C socket. If this fails, restart the device.
LED output flickers or becomes intermittent during treatment.	Loose connections	Check that the USB C connector for the face mask or neck attachment is inserted fully. Restart the device.
No sound can be heard from the controller	Volume has been muted/turned down	Select the settings icon and set the controller volume to the desired level
The display screen is too dim	Screen brightness has been set too low	Select the settings icon and set the screen brightness to the desired level
The display screen is in a different language	Wrong language has been selected	Select the settings icon and choose the desired language

If the suggested actions fail to fix the problem, or the problem reoccurs and your device is within the warranty period, please return to your provider with a proof of your purchase.

DO NOT attempt to investigate any further

The EOS-X by Déesse PRO device system contains no serviceable parts. DO NOT attempt to open and try and modify the EOS-X by Déesse PRO device in any way. Modifying the EOS-X by Déesse PRO device will void your warranty and can cause an electric shock when the system is used again.

13.0. Labelling

The labelling supplied with the device contains the following information

- Date of manufacture
- Catalogue number and batch code
- Operating environment
- Manufacturers information

14.0. Technical specification

- Model Number [REF]: TN2214
- Treatment Time: 10 minutes +/- 11 Seconds
- Wavelengths: 415nm +/- 10nm, 545nm +/- 10nm, 630nm +/- 10nm, 830 +/- 10nm and 1070m +/- 20nm

- Maximum Output Intensity: Mode 1 – Anti-Ageing (Red & Near-Infrared) 39mW/cm²
- Maximum Output Intensity: Mode 2 – Acne (Blue & Red) 53.5mW/cm²
- Operating Temperature: 0°C to 30°C
- Storage & Transport Temperature: -5°C to 55°C
- Operating Humidity: 30% to 80% without condensation
- Storage & Transport Humidity: 30% to 80% without condensation
- Electrical Input From Power Supply: 100-240Vac, 50/60Hz. 1.5A
- Power Supply Electrical Output: 24Vdc, 2A
- IP22

15.0. Explanation of Symbols

There are several safety symbols on the EOS-X by Déesse PRO device and associated packaging and in the operating manual. Use this guide to interpret them. Not following these symbols may result in one or more of the following hazards; electrical shock and/or fire.

Caution/Warning

It is mandatory that you read and comprehend all of the information in this operating manual prior to using your EOS-X by Déesse PRO device.

Refer to the operating manual for guidance on using your EOS-X by Déesse PRO Next Generation device.

The EOS-X by Déesse PRO device is specified as applied part type BF. These are parts that are applied to the patient's body, come into direct contact with it, and are designed to provide a degree of protection against electric shock.

Conforms to certain EC Directives on health, safety, and environmental protection standards. CE stands for 'Conformité Européenne'.

Conformity to UK Directives on health, safety, and environmental protection standards. UKCA stands for UK Conformity Assessed. (Wrinkles Indication Only)

Legal Manufacturers Details

Manufacturing Date

Keep Dry

Serial Number

REF	Model number for the EOS-X by Déesse PRO device: TN2214
LOT	Batch Code
MD	Medical Device
	Made in China
	Keep away from direct sunlight
	Temperature limit
	Humidity limitation
	The cardboard packaging that the EOS-X by Déesse PRO device came in can be recycled
	The outer cardboard packaging that the EOS-X by Déesse PRO device came in is FSC Certified and has been verified as being made from recycled content
IP22	The device is protected against solid objects over 12.5mm in diameter (such as fingers) and vertical water droplets when the device is tilted 150C
	Product shall not be disposed of with normal household waste (2012/19/EU)
	Input for the power supply. This means that the device may only be powered with the adapter that is supplied with the device. The usage of other adapters may not be safe and could lead to electrical shock. The adapter has the following electrical specification: 100-240Vac, 50/60Hz. 1.5A. Power Supply Electrical Output: 24V, 2A.
	Australian Sponsor (Aus only)
RP	Responsible Person (Europe only)
	This button on the user interface puts the controller into sleep mode.
	This is the ON/OFF symbol on the power switch on the rear of the controller unit (O is off, I is on)
	The device is electrically double insulated
UDI	Unique device identifier

16.0. Warranty

Under normal use, the device is guaranteed for 3 years after its original purchase date against manufacturing defects in materials or workmanship.

This warranty extends only to the original purchaser.

This warranty does not cover damage caused by the following:

- Accident, misuse, abuse or alteration
- Servicing by unauthorized persons
- Connection to incorrect power sources
- Use with unauthorized accessories

What to do if your EOS-X by Déesse PRO device is faulty

The brand owners, Harpar Grace International, or its authorised distributor shall not be responsible for any incidental, special or consequential damages resulting from the use of this product. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to one year from the date of the original purchase.

This warranty gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you.

This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model, or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose, are hereby disclaimed.

If the product is defective within the warranty period, please contact your country distributor (a list of authorised Déesse PRO distributors can be found on www.deessepro.com).

It is essential you register your warranty at www.deessepro.com/pages/warranty

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Harpar Grace International Ltd, Unit 4 Imperial Court, Magellan Close,
Walworth Business Park, Andover, Hampshire, SP10 5NT, United Kingdom

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Certainty IP Ltd, Dublin, D18 CV48. +353 1582 9003

Emergo Australia, Level 20 Tower II, Darling Park, 201 Sussex Street, Sydney,
NSW 2000. Australia

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iSmart Developments Ltd, 129 Green Lanes, Sutton Coldfield, West Midlands, B73 5LT,
United Kingdom

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