

med  the original therapeutic massager  
massager



owner's manual  
accuvibe mini

[medmassager.com](http://medmassager.com)

189 Continuum Drive, Suite A, Fletcher, NC 28732

Thank you for purchasing the Accuvibe Mini Massage Gun. In order to ensure safe and correct use of this product, please read the instructions carefully before use and operate according to instructions.

**Please consult your doctor before us if you are experiencing any of the following conditions;** Preganacy, diabetic complications (such as neorupoathy or retinal damage), wearing a pacemaker, recent surgery, epilepsy or migraine, disc herniation, spondytits, vertebral detachment or spinal joint stiffness, recent joint replacement surgery or IUD, metal or steel plate and any other serious problems related to your health.

## contents

important information .....	4
product overview .....	6
custom attachments.....	7
charging .....	10
instructions.....	10
maintenance, cleaning, storage and transport.....	11

## important information

To reduce the risk of malfunction, electric shock, injury, fire, death and/or equipment damage, please read all instructions and warnings in this manual before use.

## safety precautions

1. Minors are advised to use this product under adult supervision.
2. In case of muscle injury, please consult your doctor before using this product.
3. Do not use on wet or dirty skin. It is recommended to use over clothing.
4. Do not use over the neck, bones or body parts with less muscle coverage.  
It's only applicable on the body soft tissue without causing pain or discomfort.
5. Do not put any objects in the vent holes of this product.
6. Do not use when the body has bruises, pain or discomfort.
7. Do not immerse the product in water and keep water out of the product vents.
8. Do not put the product into a fire or close to a fire source.
9. Do not drop this product intentionally.
10. Do not disassemble this product by any means.

11. Keep hair, fingers and other parts of the body away from the massage head while in use.




**DO NOT TOUCH** this part of the product when in use.

### warnings and reminders

1. Keep away from liquids and heat sources and keep vents dust free.
2. After charging is complete, unplug the unit.
3. The product cannot be used during charging.
4. Please use this product according to the instructions.
5. Carefully check whether the product can be used normally according to the instructions in this book before each use.

## product overview

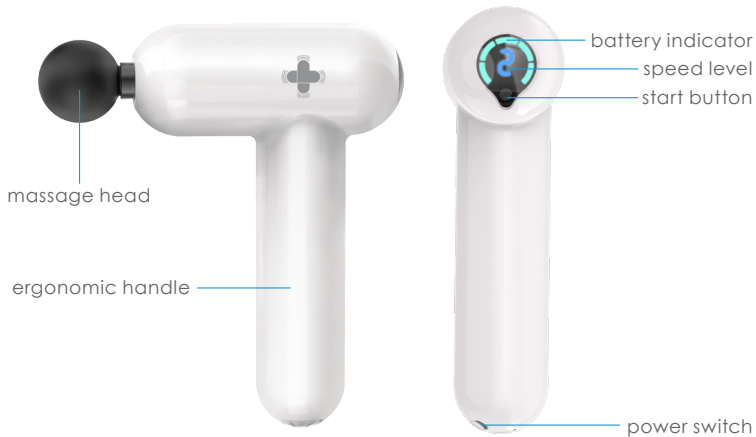

## custom massage attachments



mini achilles  
massage head



mini trigger point  
massage head



mini general  
massage head

## changing attachments

While pulling, slightly rotate the massage head to pull off. Gently push the new massage head onto the massage gun.

<b>Product Name:</b> MedMassager™ Accuvibe Mini Massager Gun	
<b>Rated Voltage:</b> 5V	<b>Rated Current:</b> about 2000mA
<b>Product Weight:</b> approx. .6kg	<b>Massage heads:</b> 3 custom attachments
<b>Power Supply:</b> Rechargeable lithium battery	<b>Mode:</b> Five speeds, manual speed control

## custom massage attachments

The MedMassager™ Accuvibe Mini Massage Gun comes with three custom massage attachments. Each massage head is specifically designed to provide a deeper, percussion massage.







### mini general massage head

The general massage head provides better fascia stimulation and is suitable for multiple muscle groups.



### mini trigger point massage head

The trigger point massage head provides deeper penetration into the trigger points to provide pain relief.



### mini achilles massage head

The achilles massage head provides massage relief for the ankles and the achilles tendon. This massage head is perfect for runners and can also be used on the back of your neck.

## charging

1. Charge the battery fully before initial use.
2. When charging, insert the USB cable into the charger plug and plug it into the wall outlet.
3. Battery level displays when charging; red flashing light (<25% battery power), yellow flashing light (<50% battery power), cyan flashing light (<75% battery power). When the light is green, the battery is fully charged.

## instructions

Turn the power switch to the ON position. Press the power button on the smart display to enter into the first speed setting. Press the power button again to enter into the next speed setting. Continue pressing the power button until the desired setting is reached. After the fifth speed level is reached, pressing the power button again will put the massage gun into stand-by mode. In any mode, hold the power button for half a second to enter into manual speed mode.

Please make sure to move the power switch to the OFF position when not in use.

## maintenance, cleaning, storage and transport

1. Clean the device with a slightly wet towel and dry with a soft cloth.
2. Place the product in the package provided while storing or transporting.
3. Store in a cool, dry, well ventilated place.
4. If not in use for a long time, store in the package after cleaning and, charge every six months.
5. This product is not suitable for use in a moving vehicle.



[medmassager.com](http://medmassager.com)

877.501.6334

1.828.475.6022 (International)

[Sales@MedMassager.com](mailto:Sales@MedMassager.com)

189 Continuum Drive, Suite A, Fletcher, NC 28732