

INSTRUCTION MANUAL

Clinician Adjustable

ASSEMBLY & OPERATION INSTRUCTIONS



Flat Top




Lift Assist Backrest Top



Lift Assist Salon Top
*Shown with optional
Adjustable Side Arm Rests
and Universal Table Extender*

made in the USA
with US & imported parts

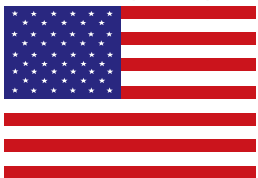

www.oakworks.com · 717.235.6807

OAKWORKS®

TABLE ASSEMBLY INSTRUCTIONS



Parts Included:



1 - Wrench




8 - Washers



8 - Nylock Nuts



4 - Main Legs w/extensions

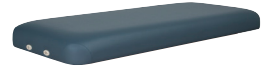

1 - Table Top



1 - Shelf

STEP 1

Place the table top upside down on the floor.

STEP 2

Place the shelf upside down on top of the table skirt. (Fig. 1) Each leg has a number by the top bolt that corresponds with a number on the inside corners of the table's base.

STEP 3

Place the top leg bolt through the metal corner bracket in one corner of the table skirt. (Fig. 1 inset) For flex top models, it may be necessary to raise the foot section to install the leg.

STEP 4

Lift the same corner of the table shelf and place the bottom leg bolt through the metal corner brackets on the shelf.

STEP 5

Place a flat washer and nut on each bolt. Do not tighten. (If cabinet is included, please go to step 1 below.)

STEP 6

Repeat steps 3-5 on remaining corners, installing the remaining three leg assemblies to the skirt and shelf.

STEP 7

Using the included wrench, or a 1/2" wrench, tighten all eight nuts securely. When tight, assembled unit should look like figure 2.

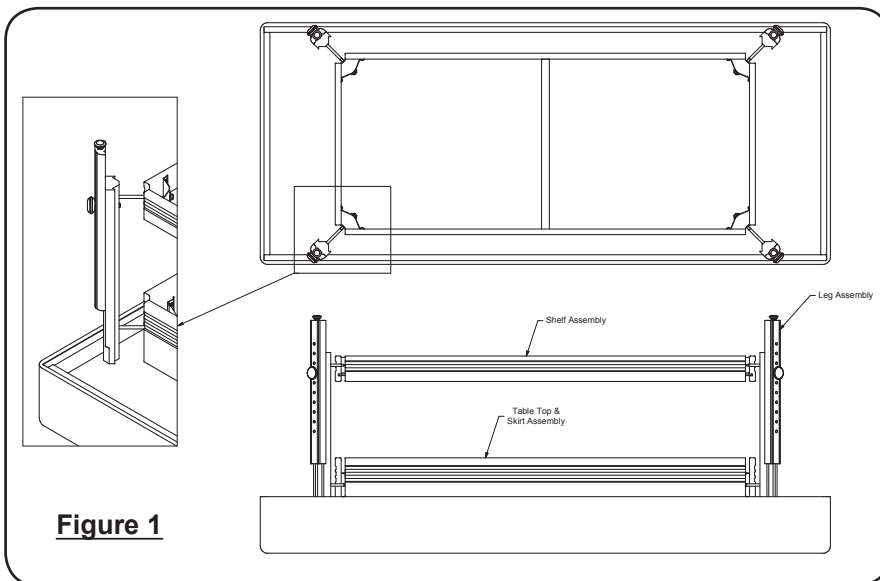


Figure 1


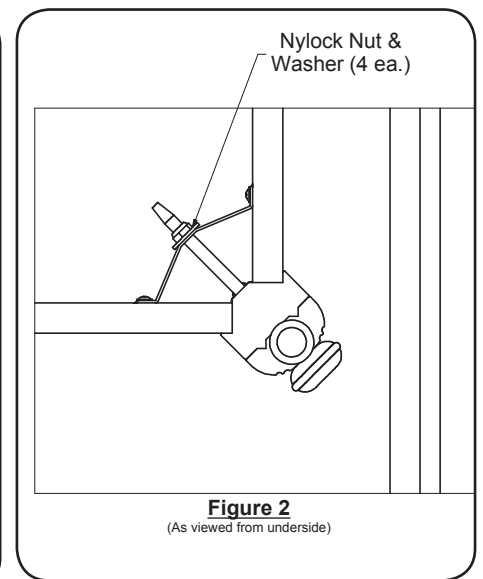


Figure 2
(As viewed from underside)

CABINET OPTION INSTALLATION INSTRUCTIONS

STEP 1

Tighten all four nuts on the bottom shelf, and upright the table with assistance from another person.

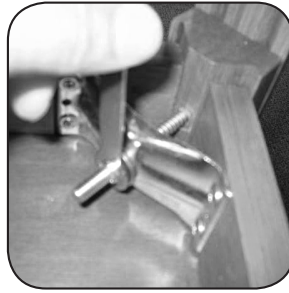


STEP 2

Place end panel in position. Insert supplied screws (2) into the hole in the legs on both sides of the end panel. Tighten screw until end panel is secured. Do not over tighten. Tighten just enough to hold end panel in position. Repeat for other end panel.




STEP 3

Remove lower door track on shelf with supplied screw driver. Lay doors on shelf in the position they will be installed.

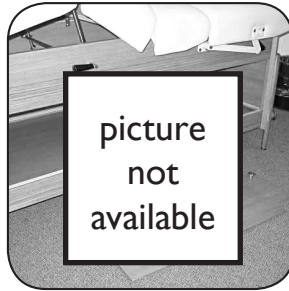

STEP 4

At an angle, install the top of the two doors into the upper door track. Place the lower track on the bottom of the two doors. Slide the lower track into the final location and replace the screws in the lower track.

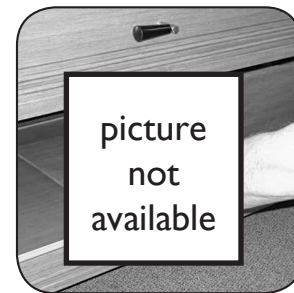

STEP 5


Slide the doors back and forth to insure easy movement.

Repeat on other side.



OPERATING THE MANUAL TOP (IF EQUIPPED)

BACK REST

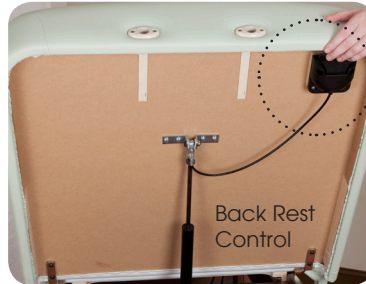

To Raise the Back Rest:

Raise the backrest by pressing the lever under the head section and lifting up to the desired height. It is recommended that the client takes some of their weight off the backrest while you lift.

To Lower the Back Rest:

Press the lever under the head section and push the backrest down. Release the lever and the backrest will stay in place.

To change position of the backrest lever:

The backrest lever can be moved to the opposite side of the backrest by removing the screws and relocating the hand lever to the opposite side of the top into the four pre-drilled holes.

LEG REST



To Raise the Leg Rest:

Raise the leg rest by pulling the handle on the side of the leg section and raise to the desired height. It is recommended that the client takes some of their weight off the leg rest while you lift.

To Lower the Leg Rest:

Pull the handle on the side of the leg section and push the leg section down. Release the handle once leg rest is at desired location and the leg section will stay in place.

To change position of the leg rest handle:

The leg rest handle can be moved to the opposite side of the table. First, remove the screws from the handle and replace just the screws. Move handle to the opposite side of the table. Remove the screws that are there and mount handle.

ARM REST INSTALLATION AND ADJUSTMENT INSTRUCTIONS

OPTIONAL ACCESSORY AVAILABLE ON SALON TOP ONLY

CAUTION: Do not use the arm rest for customer weight support during the mounting or dismounting of the table. The arm rest is not designed to support the patient's weight.



Adjustable Side Arm Rests



Manicure Arm Rests

TO INSTALL THE ARM REST:

The arm rests adjust automatically as the table's back support moves to remain parallel to the floor in any position. They are designed to be strong and will take a reasonable amount of weight, but they are not designed to be used to sustain a person's weight when getting on or off the table.



1. Insert end of arm rest into hole in backrest while holding arm rest at 60° angle.



2. Push end all the way into hole and drop arm rest into horizontal position. Tug slightly to make sure the armrest is fully locked to the table (for patient safety).



To remove the arm rest, simply reverse the procedure above.

TO LEVEL THE ARM REST:


The arm rest height is adjusted to level at the factory. When in use, some of the new parts will wear, and you may see some sagging of the arm rest over time. If the arm rest is too low the height can easily be adjusted using the enclosed Allen wrench.



To Raise the Arm Rest



1. Loosen bottom 2 set screws (turn counter-clockwise).

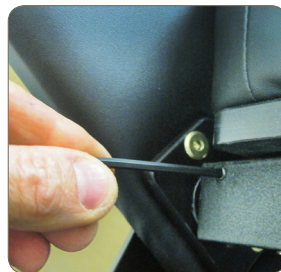

2. Turn side set screw in (clockwise) while lightly raising the arm rest with your other hand (to make the adjustment easier).



3. When you have achieved the desired height, tighten the two bottom set screws (turn clockwise).

To Lower the Arm Rest

1. Loosen bottom 2 set screws (turn counter-clockwise).

2. Turn side set screw out (turn counter-clockwise).

3. When you have achieved the desired height, tighten the two bottom set screws (turn clockwise).

THIS PAGE INTENTIONALLY LEFT BLANK

NOTES:

THIS PAGE INTENTIONALLY LEFT BLANK

NOTES:

Clinician Adjustable

ASSEMBLY & OPERATION INSTRUCTIONS

CONTACT INFORMATION:

OAKWORKS® Inc.

923 East Wellspring Road
New Freedom, PA 17349

Phone: 717-235-6807

FAX: 717-235-6798

www.oakworks.com



European Compliance for CE Information:
Emergo Europe
Molenstraat 15
The Hague, 2513 BH Netherlands
info@emergogroup.com

CAUTION

Maximum Working
Load Rating:
500 LBS. (227 KG.)

Manual Part Number MMMNST0023


Revision date: 02.01.12

Edition 1 - February 2012

Printed in English

Printed in USA

made in the USA
with US & imported parts

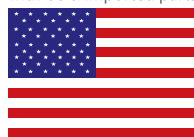

OAKWORKS®