

MODELS 4418 AND 4416 SPA PRO™ MASSAGE STONE HEATER

Introduction to Stone Message Therapy

Stone massage therapy is a wonderful procedure that can bring many therapeutic benefits when properly administered. As with any therapy, proper care must be practiced to ensure the safety of the individual receiving the therapy. Before using the SPA PRO™, we strongly recommend instruction on stone massage therapy from a qualified, professional therapist.

Table of Contents

IMPORTANT SAFEGUARDS	3
SAVE THESE INSTRUCTIONS	3
CAUTIONS	4
BEFORE USING YOUR SPA PRO™ FOR THE FIRST TIME	4
DIRECTIONS FOR USE	5
MAINTENANCE AND CLEANING	6
GET TO KNOW YOUR SPA PRO™ MASSAGE STONE HEATER	6-7
LIMITED WARRANTY AND SERVICE SATISFACTION	8
FIVE EASY STEPS TO SATISFACTION	9



The Metal Ware Corporation
1700 Monroe St., P.O. Box 237
Two Rivers, WI 54241-0237
Phone: (800) 288-4545
FAX: (920) 794-3164
Website: www.nesco.com
Part Number: 63943

©The Metal Ware Corporation., 2003

IMPORTANT SAFEGUARDS

When using the SPA PRO™ or any electrical appliances, these basic safety precautions should be followed, including the following:

1. **READ ALL INSTRUCTIONS BEFORE USING APPLIANCE.**
2. Do not touch hot surfaces. Use handles or knobs.
3. To protect against electrical hazard, do not immerse cord, plug or base unit in water or other liquid.
4. Always plug appliance cord into the wall outlet, then adjust Temperature Control to minimum (off), then remove plug from wall outlet.
5. Do not operate when bathing or in the shower.
6. Close supervision is necessary when any appliance is used by or near children.
7. Unplug from the outlet when not in use and before cleaning. Allow to cool before putting on or taking off parts, and before cleaning the appliance.
8. Do not operate any appliance with a damaged cord or plug or after the appliance malfunctions or has been dropped or damaged in any manner. Return appliance to The Metal Ware Factory Service Center for examination and repair or electrical/mechanical adjustments.
9. The use of accessories or attachments not recommended by the manufacturer may cause fire, electrical shock or injury.
10. Do not use outdoors.
11. Do not let cord hang over the edge of the table or counter, or touch hot surfaces.
12. Do not place on or near a hot gas or electric burner or in a heated oven.
13. Extreme caution must be used when moving an appliance containing hot liquids.
14. Do not use appliance for other than intended use as described in this manual.
15. Do not fill the water reservoir to more than 1 inch below the top rim. Place appliance on a level surface.
16. Always place the stones in the water reservoir first, then add water.
17. Do not use if you have rashes, burns, cuts or other skin irritations, or if you have sensitive or damaged skin.
18. Never place an appliance in a tub or sink to fill while it is plugged into a wall outlet.
ELECTROCUTION HAZARD!
19. Always test the temperature of each stone by touching it with your fingers first and hold it in your hand before placing the heated stone on another person.
20. Use heat resistant gloves or slotted wooden spoon to lift heated stones from the SPA PRO™ to avoid burns or scalding.
21. Never add boiling or steaming water to avoid scalding and/or burns.

SAVE THESE INSTRUCTIONS

CAUTION: A short power supply cord is provided to reduce the risk resulting from becoming entangled in or tripping over a long cord.

CAUTIONS

1. The SPA PRO™ Massage Stone Heater should not be used if the following conditions are present: Tuberculosis, tumors, open wounds, hemorrhages, bruises, broken skin, thrombosis or phlebitis, varicose veins, or undetermined joint pain or pain in the calf.
2. People who are pregnant, have diabetes or other illness, should have a physician's approval prior to receiving therapy.
3. Contact your physician if you have any questions concerning the therapeutic affects of massage stones heated with this appliance.
4. Discontinue use of the SPA PRO™ if you experience any discomfort from the therapy and contact your physician.
5. Do not use this appliance when you are feeling drowsy or are going to fall asleep.
6. The SPA PRO™ is not a toy. **USE EXTREME CAUTION** when using this appliance and **PLEASE** keep children away from this appliance.

Before Using Your SPA PRO™ For The First Time

1. Unpack your SPA PRO™ Massage Stone Heater completely and REMOVE THE WATER RESERVIOR, ALL PACKAGING MATERIAL AND INFORMATION LABELS. Place unit in a WELL VENTILATED AREA ON A HEAT RESISTANT SURFACE.
2. Turn the Temperature Control full counter-clockwise and plug the SPA PRO™ into a suitable 120 volt AC outlet.
3. Turn the Temperature Control to the "Preheat" setting and allow the SPA PRO™ to operate for approximately 60 minutes, or until any odor or smoking disappear.
4. Some smoke and odor are a normal part of curing the heating element and will not occur again during normal spa use.
5. When curing is completed, turn the Temperature Control full counter-clockwise to turn the appliance off, unplug SPA PRO™ from wall outlet and allow unit to cool.

DIRECTIONS FOR USE


1. Wash all stones before each use.
2. Place a hand towel on the bottom of the SPA PRO™ water reservoir. Place all stones on top of the towel.
3. Before plugging appliance into outlet, add sufficient water fill the water reservoir to just cover stones or fill to maximum capacity, but no more than 1-inch below top rim. Turn the Temperature Control full counter clockwise. Place Cover on appliance.
4. Plug power cord into a standard 120 VAC grounded outlet. Turn the Temperature Control to "Preheat" position for about 20 to 30 minutes to speed warm up time. When the water reaches the desired temperature, turn the Temperature Control to the "Operating Range" to maintain the proper temperature.
5. Use an accurate thermometer to measure water temperature in the Water Reservoir to ensure that stones are at a proper therapeutic temperature. Consult a professional therapist for correct stone temperatures if you are unsure.
6. Model 4418 (18 qt version) only: The Indicator Lamp cycles off and on with the heating element. The element is heating when the lamp is on.
7. Carefully remove the first stone from the SPA PRO™ Water Reservoir. Test the temperature of the stone by touching it with your finger, then hold it in your hand. The temperature should feel warm to the touch, but not uncomfortable.
8. When the therapy session is finished, turn the Temperature Control full counter clockwise to turn unit off. Always unplug appliance when not in use.
9. Do not drop stones into the SPA PRO™ heating unit. This may damage the appliance. Use care when removing or replacing stones.
10. Do not apply warm stones to areas of your body that are inflamed.

Maintenance and Cleaning


1. Clean unit after each therapy session.
2. Turn Temperature Control to minimum and unplug SPA PRO™ from wall outlet before cleaning.
3. Allow appliance to cool to room temperature before cleaning or add or removing parts.
4. Carefully remove all stones first, then remove the water.
5. Clean Water Reservoir and Cover with a solution of warm soapy water, then rinse with clean water and dry with a towel.
6. Wipe Main Heating Unit with a damp cloth. Special care must taken when cleaning the Main Heating Unit. **NEVER** immerse Main Heating Unit in water or place in a sink to rinse.
7. Rinse the stones in warm water with a small amount of mild detergent added. Rinse stones with clean warm water and allow to dry before storing.

Get To Know Your SPA PRO™ Massage Stone Heater

Model 4416 – 6 qt. SPA PRO™

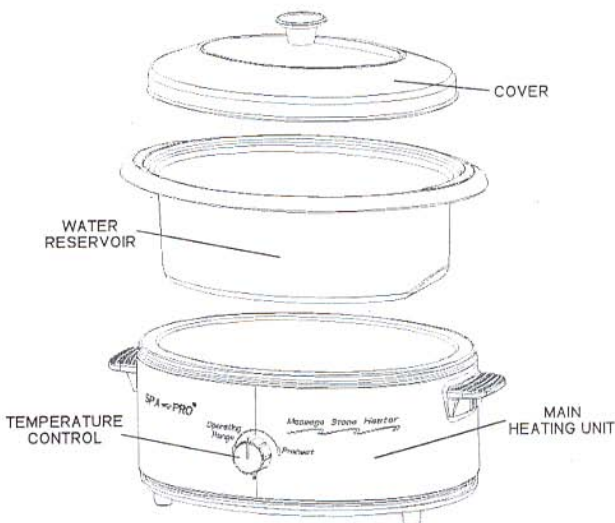

Model 4418 - 18 qt. SPA PRO™

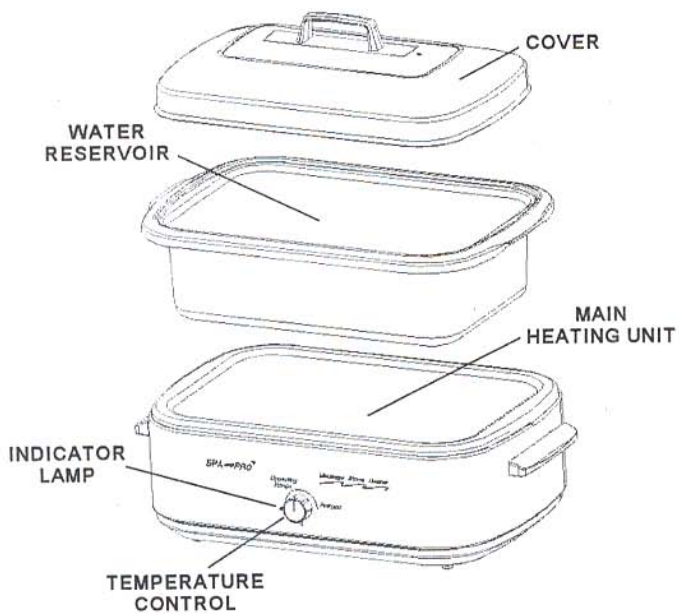

ONE YEAR LIMITED WARRANTY

This appliance is warranted for one year from date of purchase against defects in material and workmanship. This warranty does not cover transportation damage, misuse, accident or similar incident. This warranty gives you specific legal rights and you may have other rights, which vary from state to state.

Your new SPA PRO™ Massage Stone Heater comes equipped with numerous safety features. Any attempt to interfere with the operation of these safety features makes this warranty null and void. In the event we receive a SPA PRO™ Massage Stone Heater for service that has been tampered with, we reserve the right to restore it to its original state and charge for the repair.

Defective products may be returned, postage prepaid, with a description of the defect to: The Metal Ware Corporation, 1700 Monroe Street, Two Rivers, Wisconsin 54241 for no charge repair or replacement at our option.

Service and genuine replacement parts may be obtained from The Metal Ware Corp. Factory Service Department.

For service in warranty, follow instructions set forth in warranty. When ordering new parts, make sure that you always mention the model number of the product, which is found on the back of the unit.

About Your Warranty and Service

One simple step sets you on the road to complete satisfaction:

Mail in your warranty card and file your SPA PRO™ Massage Stone Heater's exclusive registration number now!

Your serial number and SPA PRO™ Massage Stone Heater Model are located on the label at the back of the appliance.

We enter this number into our database. It registers your warranty and identifies your SPA PRO™ Massage Stone Heater specifically by its individual code. Should you call our 800 toll free line, your registration number helps us identify your appliance and speeds up our assistance to you.

Five Easy Steps to Satisfaction

Should you have a problem with your SPA PRO™, please refer to these steps for assistance:

1. Call us on our toll free number, 1-800-288-4545 and tell use about your problem.
2. If we instruct you to send all or part of your appliance to us for repair or replacement, the Customer Satisfaction representative will give an authorization number (this number notifies our Receiving Department to expedite your repair). Pack your unit carefully in a sturdy carton to prevent damage because any damage caused by shipping is not covered by the warranty.
3. Print your name, address and authorization number on the carton.
4. Write a letter explaining the problem. Include the following: your name, address and telephone number and a copy of the original bill of sale.
5. Attach the sealed envelope containing the letter inside the carton. Insure the package for the value of the SPA PRO™ Massage Stone Heater and ship prepaid to:

**The Metal Ware Corporation
1700 Monroe Street
Two Rivers, WI 54241**