

# INSTRUCTION MANUAL

---

**OAKWORKS®**

---

## *Relaxation Lounge™*



Manual Part # MMMNST0007  
Edition 1, English Sept 2007  
Edition 2, English May 2008

All rights are reserved. No part  
of this document may be  
photocopied, reproduced or  
translated to another language  
without prior written consent  
of OAKWORKS®, Inc.

OAKWORKS® is a registered  
trademark of OAKWORKS®, Inc.

#### Notice

The information contained in this  
document is subject to change  
without notice and should not  
be construed as a commitment by  
OAKWORKS®, Inc.

OAKWORKS®, Inc. assumes no  
responsibility for any errors that  
may appear in this document nor  
does it make expressed or implied  
warranty of any kind with regard  
to this material, including, but not  
limited to, the implied warranties  
of merchantability and fitness for  
a particular purpose.

OAKWORKS®, Inc. shall not be liable  
for incidental or consequential  
damages in connection with or  
arising out of the furnishing,  
performance, or use of this  
document and the program  
material which it describes.

---

## TABLE OF CONTENTS

Symbol Identification .....	pg ii
Warnings .....	pg iii-v
Grounding Instructions .....	pg v
Introduction .....	pg 1
<b>Section I: Product Description &amp; Photo</b>	
<i>Relaxation Lounge</i> <sup>™</sup> .....	pg 2-3
<b>Section II: Installation</b>	
Leveling Instructions .....	pg 4
<b>Section III: Directions for Use</b>	
Moving the <i>Relaxation Lounge</i> <sup>™</sup> .....	pg 5
Installing the Side Console .....	pg 6
Installing the optional reading light .....	pg 9
Using the Optional Sound System .....	pg 10
<b>Section IV: Troubleshooting</b>	
If Table Will Not Change Position .....	pg 11
<b>Section V: Care &amp; Cleaning</b> .....	pg 12
<b>Section VI: Product Identification</b> .....	pg 13
<b>Section VII: Table Specifications</b> .....	pg 14
<b>Section VIII: Warranty</b> .....	pg 15
<b>Section IX: List of Parts</b> .....	pg 16
<b>Contact Information</b> .....	back cover

## SYMBOL IDENTIFICATION

### SYMBOL IDENTIFICATION



This symbol, when used in documentation on this product, represents a caution. Be sure to read and comply with all precautions and warnings.



This symbol, when used in documentation on this product, indicates danger and warns against an electrical shock hazard. Be sure to observe and comply with all warnings and instructions.



This symbol, when used in documentation on this product, indicates a Protective Earth (Ground) Terminal.



This symbol when used in documentation on this product, warns that during transport there should be no stacking of containers.



This symbol, when used in documentation on this product, indicates that the product should be protected from moisture. The humidity specifications for Transport & Storage are listed on page 14.



This symbol, when used in documentation on this product indicates that information is given regarding the recommended temperature limits during transport and storing. The specifications for these limits are listed on page 14.



This symbol, when used in documentation on this product, indicates the date of manufacture of the device.



This symbol, when used in documentation on this product indicates alternating current (AC).



This symbol, when used in documentation on this product, indicates direct current (DC).



This symbol, when used in documentation on this product, indicates that instructional information is included. All instructions should be read carefully and completely before using the equipment.

## WARNINGS & IMPORTANT SAFETY INSTRUCTIONS

READ AND  
SAVE THESE  
INSTRUCTIONS

### IMPORTANT SAFETY INSTRUCTIONS (CONT.)

#### WARNING

Close supervision is necessary when this furnishing is used by, or near, children, invalids or disabled persons. DO NOT leave children alone or allow them to operate the controls. Keep controls out of reach of children. DO NOT allow children to crawl into the cabinet or under the table.

#### WARNING

The table weighs as much as 250 lbs. DO NOT lift or move the unit without using proper safety procedures or equipment. DO NOT turn the *Relaxation Lounge*™ on its side for any reason. To do so could cause damage to the base and will void the warranty. There are no serviceable items that require this action. Never lift the table by grasping the table top as the table top is not designed to hold the weight of the base and you could damage the table.

#### CAUTION

When operating the *Relaxation Lounge*™ controls and making any and all adjustments to the table, move hands and feet away from the working mechanism of the table.

#### CAUTION

DO NOT store or place anything on the cabinet hood under the backrest as this can cause damage to top and/or hood and will void the warranty.

#### CAUTION

Maximum Working Load Rating (client and accessories): 500 lbs.

## WARNINGS & IMPORTANT SAFETY INSTRUCTIONS

READ AND  
SAVE THESE  
INSTRUCTIONS

### IMPORTANT SAFETY INSTRUCTIONS



Improper use of this device can cause injury. Be sure to read all operating instructions prior to use to prevent injury to operator or client.

DANGER



*Electrical Shock Hazard.* To reduce the risk of electric shock, always disconnect the table from the electrical outlet before moving the table or doing any maintenance or troubleshooting procedures. Never use this equipment if it has a damaged cord or plug. Use this equipment only for its intended use as described in these instructions. Do not use any attachment or accessories not recommended by the manufacturer.

WARNING



To reduce the risk of burns, fire, electrical shock or injury:

- The power supply/control module is located under the top. Refer servicing to qualified personnel and follow manual instructions when doing any authorized troubleshooting. Before doing any work on the unit, disconnect the power cord from the power outlet.
- Before cleaning with any liquid cleaner be sure to un-plug the power cord from the power outlet. Do not operate the table in standing water.
- Before moving the table, disconnect the power cord from the power outlet and secure the electrical cord, out of the way for moving.
- *Risk of Electrical Shock.* — Connect this table to a properly grounded outlet only. See Grounding Instructions.
- Use this furnishing only for its intended use as described in these instructions. Do not use attachments not recommended by the manufacturer.
- Never operate this table if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the table to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

WARNING



Risk of Injury — Keep children from moving headrest and foot support.

WARNING



Never drop or insert any object into any opening.

## WARNINGS & IMPORTANT SAFETY INSTRUCTIONS

READ AND  
SAVE THESE  
INSTRUCTIONS


### GROUNDING INSTRUCTIONS

#### IMPORTANT SAFETY INSTRUCTIONS (CONT.)

##### GROUNDING INSTRUCTIONS:

This product must be grounded. If it should malfunction or breakdown, grounding provides a path of least resistance for electrical current to reduce the risk of electric shock. This product is equipped with a cord having an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances (see example below of a nominal 120-volt circuit). Make sure that the product is connected to an outlet having the same configuration as the plug.

*No adapter should be used with this product.*

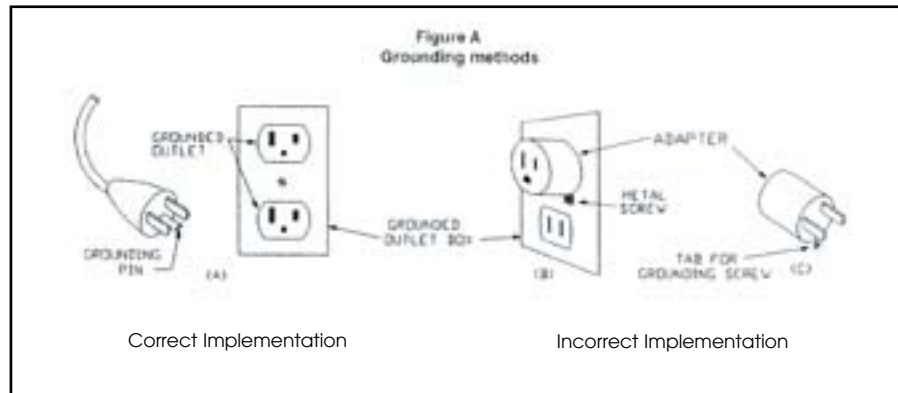

**DANGER**-Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or service man if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product - if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

**WARNING**

The *Relaxation Lounge*™ should not be used in any other manner than that described within this manual. Use of any unnecessary force or unauthorized procedures could result in injury to practitioner, client or equipment.

---

## INTRODUCTION

# *Relaxation Lounge*™

### 500 lbs. Maximum Working Load

Electronically adjustable top, controlled by a multifunctional hand control, allows you to effortlessly change your position while relaxing on the lounge.


The *Relaxation Lounge*™ should not be used in any other manner than that described within this manual. Use of any unnecessary force or unauthorized procedures could result in injury to practitioner, client or equipment.

The variety of reclining positions align the *Relaxation Lounge*™ top so your body is at the optimum level for comfort. It features a single hand control to make the precise adjustments to the back and knee rest positions allowing a custom fit to each client.

## PRODUCT DESCRIPTION

### *Relaxation Lounge™*

Model REMPST

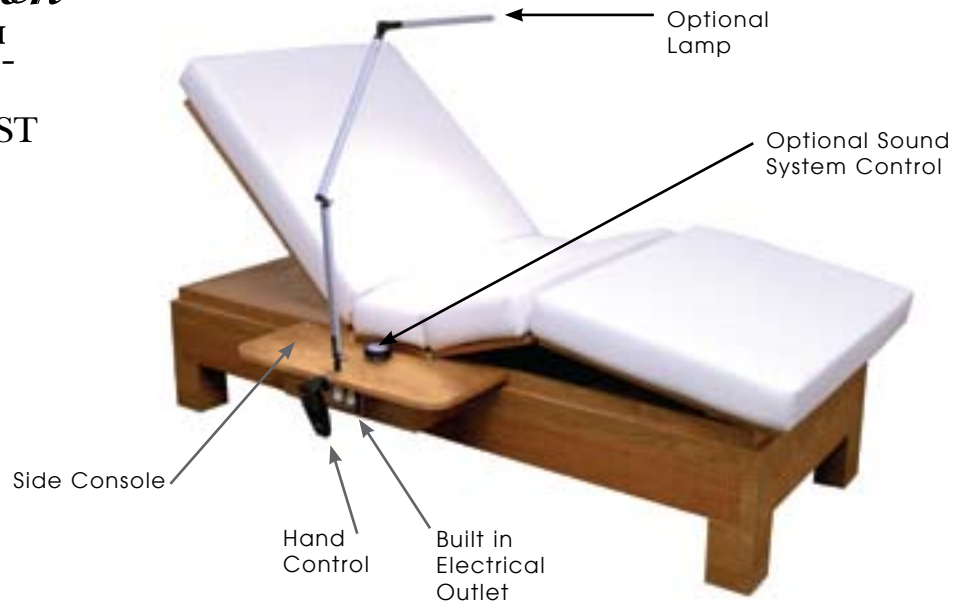

#### TABLE SPECS\*:

- Height \*: 19.5"
- Width: 32"
- Length: 78"
- Weight: 250 lbs.
- Top: Electronically Controlled Flex Top (Backrest angle 0-55°, leg rest height 0"-9")
- Electronics: 1 Hand Control
- Padding: Aero·Cel™ System
- Upholstery: Ultratouch II or TerraTouch™
- Cabinetry: Wood/Maple and Birch

For electrical specs, see Specifications Chart, pg 14.

\*Measurements are determined from floor to top of seat pad using standard padding; these measurements can be affected by natural material variations in padding

\*\*Published weight may differ due to options chosen, natural material variations and padding density

## DIRECTIONS FOR USE

- Adjusting Height and Top Positions of the *Relaxation Lounge™* using the Hand Control

## DIRECTIONS FOR USE

### *Relaxation Lounge™*

#### • Hand Control



When operating the *Relaxation Lounge™* controls be sure to observe all cautions and warnings given in the manual to prevent injury to both operator and client.




### CAUTIONS

Do not sit on the elevated Back Rest or Leg Rest.

Children must be supervised when around this equipment by a responsible adult. Do not leave children alone or allow them to operate the controls. Do not allow children to crawl into the cabinet or under the table.

### ▶▶ IMPORTANT

## HAND CONTROL OPERATIONS • *Relaxation Lounge™*

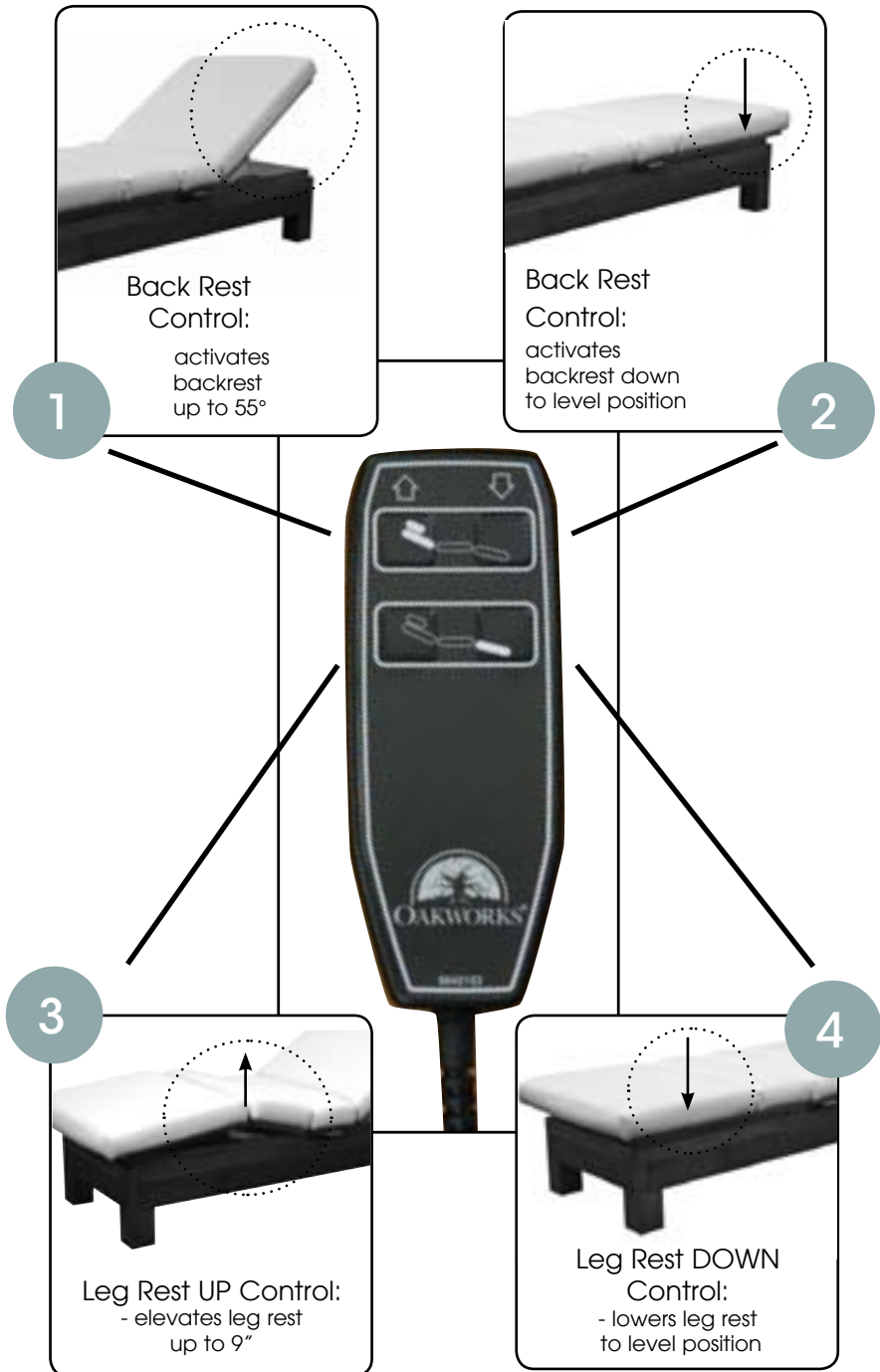

Press and hold the appropriate button until the position desired is achieved and then release.

## INSTALLATION

### • Leveling Feet



When operating the *Relaxation Lounge™* controls be sure to observe all cautions and warnings given in the manual to prevent injury to both operator and client.

The *Relaxation Lounge™* comes completely assembled and ready to use. Plug the cord into a functioning outlet that is rated for the table.

Arrange the power cord and control cords so that they will not create a tripping hazard. Place the controls in a safe location where they are convenient and accessible. Keep the controls out of the way of all table mechanisms.

### LEVELING INSTRUCTIONS:

#### If Lounge Will Not Sit Flat on the Floor:

The *Relaxation Lounge™* is equipped with Leveler Feet that adjust to level the table on uneven floors. To adjust the Leveler Feet to suit common floor differences you will need to use a 1/2" open-end wrench. Adjust each leg individually in turn by adjusting the height of each leveler foot until the unit rests level and steady. This adjustment will not effect the function of the table.



### ▶▶ IMPORTANT



### CAUTION

The table weighs as much as 250 lbs. **DO NOT** lift or move the unit without using proper safety procedures or equipment. **DO NOT** turn the *Relaxation Lounge™* on its side for any reason. To do so could cause damage to the base and *will* void the warranty. There are no serviceable items that require this action.

Never lift the table by grasping the table top as the table top is not designed to hold the weight of the base and you could damage the table.

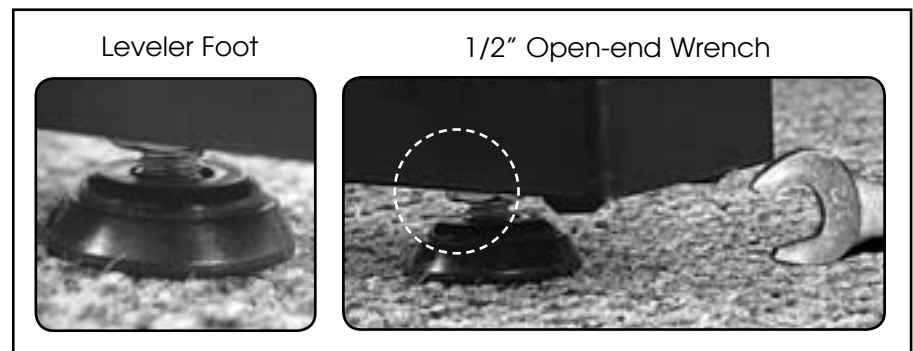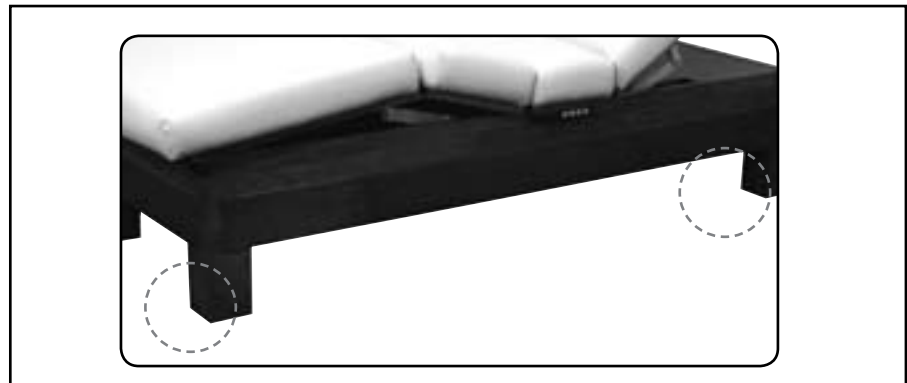

## DIRECTIONS FOR USE

- Moving the *Relaxation Lounge™*

## DIRECTIONS FOR USE

- Moving the *Relaxation Lounge™*



### CAUTION


The lounge weighs 250 lbs. **DO NOT** lift or move the unit without using proper safety procedures or equipment. **DO NOT** turn the *Relaxation Lounge™* on its side for any reason. To do so could cause damage to the base and *will* void the warranty. There are no serviceable items that require this action. Never lift the lounge by grasping the lounge top as the table top is not designed to hold the weight of the base and you could damage the table.



### CAUTION

Before moving the *Relaxation Lounge™*, disconnect the power cord from the power outlet and secure the electrical cord, hand control and foot controls for moving.

## MOVING THE *Relaxation Lounge™*

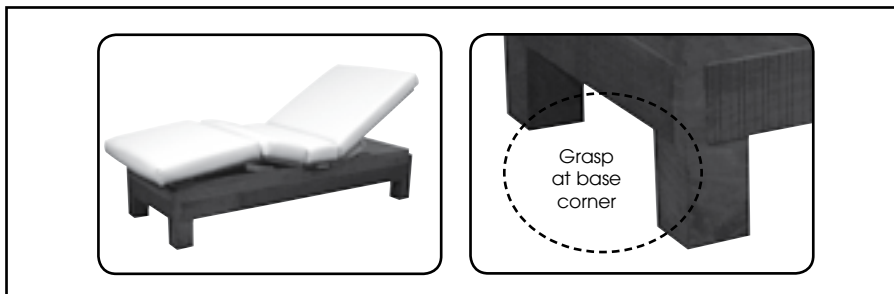

The *Relaxation Lounge™* weighs approximately 250 lbs. You will need 4 people working together if it is necessary to move the table. Never lift the table by grasping the table top as the table top is not designed to hold the weight of the base and you could damage the table and the electronics that run the tabletop.

If you need to move the table, place one person at each corner of the table. We recommend that you have the table top in a flex position to have the table top ends out of the way as much as possible. Have each person firmly grasp the cabinet base at the bottom by the table legs and lift the table from the base. As when moving any piece of heavy equipment, bend at the knee and lift with your legs beneath you and not with your back. Be aware and use proper, safe procedures for moving a table of this size to prevent injury to yourself or the table. We strongly recommend that a flatbed dolly be used that can hold the entire table if moving the table involves any distance.

## ▶▶ IMPORTANT




### CAUTION

**DO NOT** store or place anything on the hood under the backrest as this can cause damage to the top and/or hood.

## IMPORTANT SAFETY NOTE • *Relaxation Lounge™*

The hood area under the backrest top is not designed for storage or as an area to hold accessories while with a client. **DO NOT** use this area for holding *anything*. Damage can occur to the table, especially when commands that are accessed by the hand control are activated. Placing anything on this area will prohibit movement of the backrest and can cause damage to the table or the object and will void the warranty.

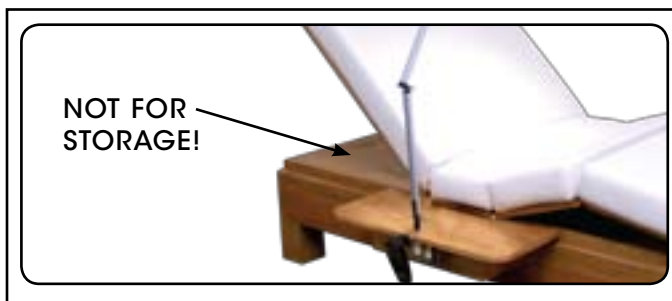


## DIRECTIONS FOR USE


## INSTALLING THE SIDE CONSOLE

Materials Needed: 1 shelf assembly  
1 hanger assembly  
1 lock chip


1. Remove the pads from the Relaxation Lounge.
2. Install the hanger assembly as shown, straddling the seat board on either side of the lounge.

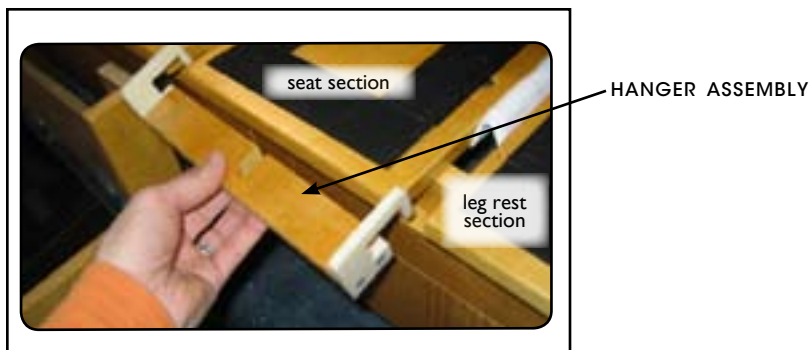

3. Push hanger all the way down as shown.

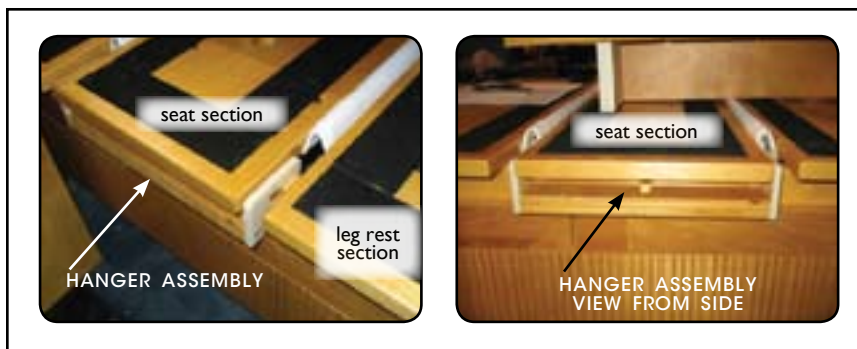

4. Insert shelf between the seat board and the hanger assembly as shown. Align notched edge with hanger assembly. Push all the way back against the lounge.

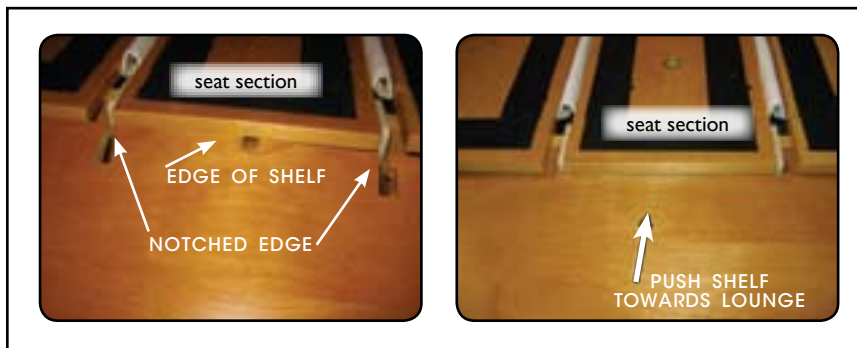


## DIRECTIONS FOR USE

- Installing the side console for the *Relaxation Lounge™*

## DIRECTIONS FOR USE

### INSTALLING THE SIDE CONSOLE (CONTINUED)

5. Push the shelf to the left to lock on the hanger assembly.

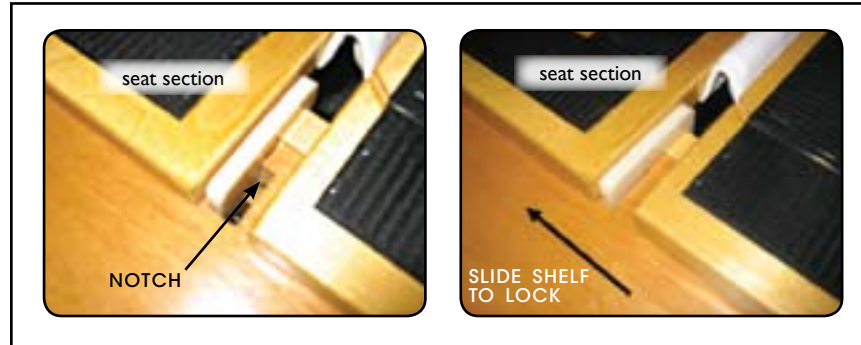

6. Insert Lock Chip in the gap on the left side as shown.




7. Feed the wires under the relaxation lounge and clip them into the provided wire clips.




## DIRECTIONS FOR USE

8. Remove the backrest dust cover and lift out the wire grommet cap. Place the side console power cord through the grommet system and plug the power cord into the provided power strip inside the table.

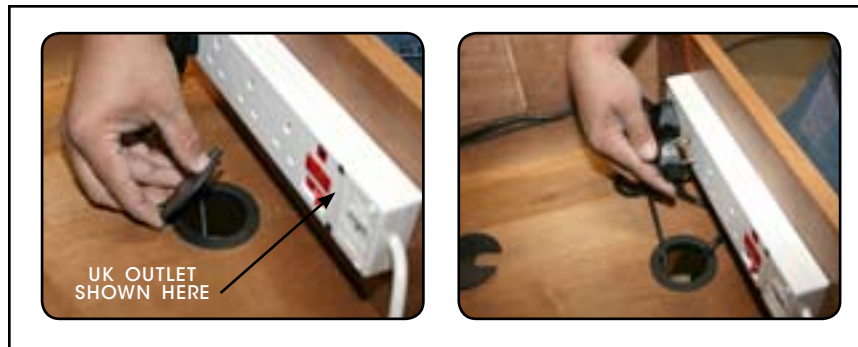

9. Replace the grommet cap, dust cover and Relaxation Lounge pads.

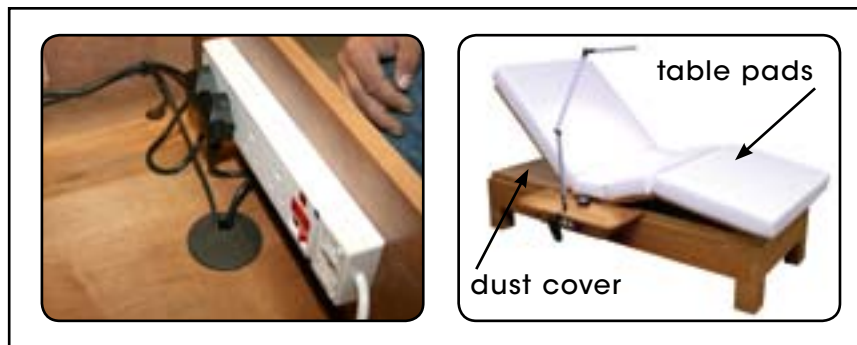

## DIRECTIONS FOR USE


- Installing the optional reading light

## DIRECTIONS FOR USE


- Optional Reading Light

### INSTALLING THE OPTIONAL READING LIGHT


1. Remove the black cap from the shelf. Insert the pin provided in the light kit so that it is flush with the under side of the shelf.

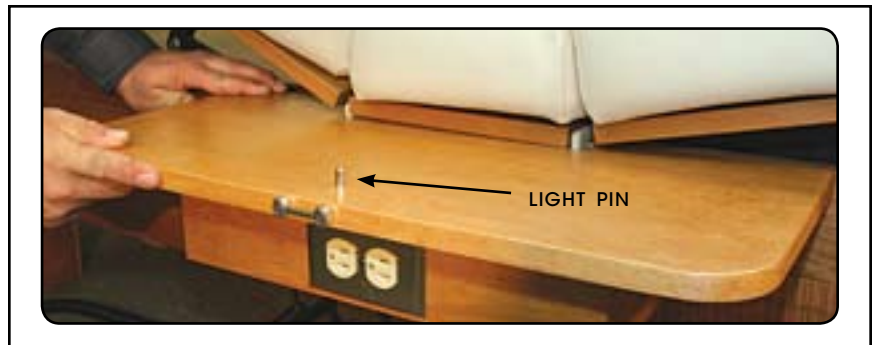

2. Insert the screw in the pre-drilled hole on the under-side of the shelf.

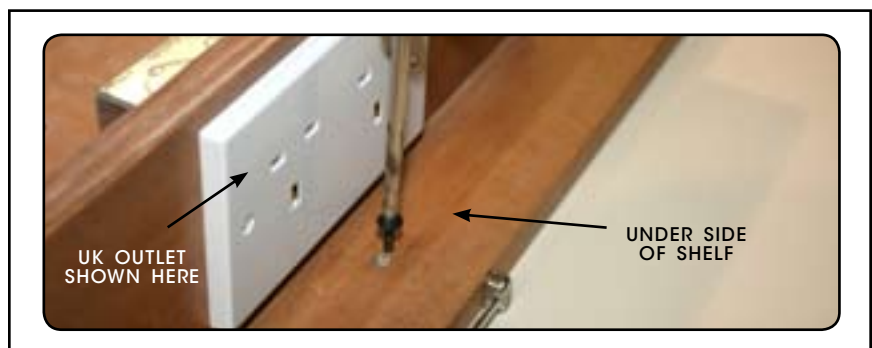

3. Remove the reading light from the box and place into the pin in the shelf as shown below.

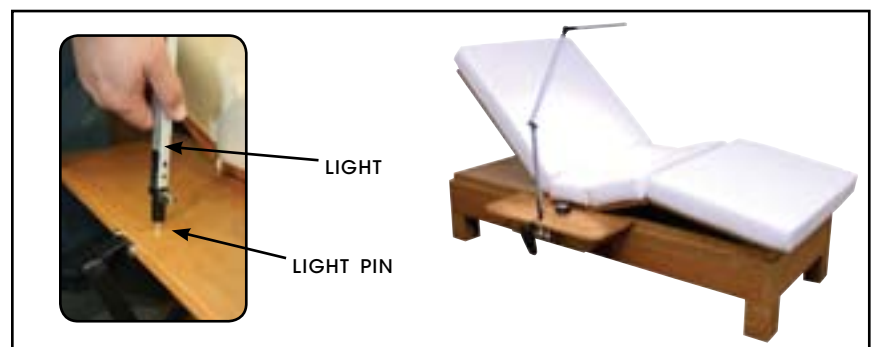


4. Plug the power cord into the outlet.

## DIRECTIONS FOR USE

- Optional Sound System

### OPTIONAL SOUND SYSTEM

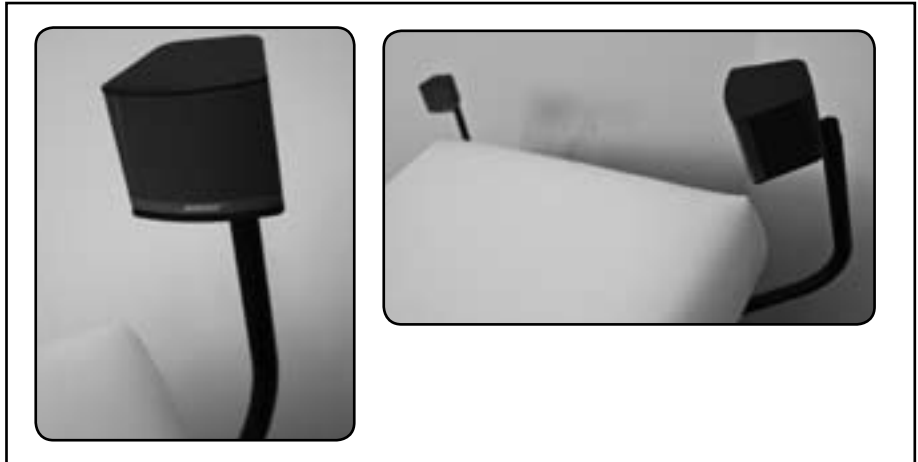

A personal sound system is available for the *Relaxation Lounge*™. The system consists of a 3 speaker system with 2 speakers located at the head end of the lounge, and one speaker located inside the lounge or under the optional console. The speakers on the head end of the lounge can be moved out of the way when not in use. The system can play music from MP3 or any source with a headphone output. The volume control and on/off switch is conveniently located next to the table control buttons.

The MP3 input can be plugged into the volume control or direct into the base speaker unit (Sometimes located inside the frame) Access is obtained by removing the headrest dust cover.



### ►► IMPORTANT



### CAUTION

Note: To remove dust cover, slide forward 1/2 inch, then lift.

Caution - When the speakers are folded in the down position make sure they are clear of the frame before operating any controls! (Speakers in contact with the frame can be severely damaged!)

## TROUBLESHOOTING

- If *Relaxation Lounge™* will not Change Height or Position

## TROUBLESHOOTING

- Controls

### ▶ IMPORTANT



### CAUTION

The lounge weighs 250 lbs. **DO NOT** lift or move the unit without using proper safety procedures or equipment. **DO NOT** turn the *Relaxation Lounge™* on its side for any reason. To do so could cause damage to the base and *will* void the warranty. There are no serviceable items that require this action. Never lift the lounge by grasping the lounge top as the table top is not designed to hold the weight of the base and you could damage the table.

### IF THE *Relaxation Lounge™* WILL NOT CHANGE POSITION:

- Check the outlet to be sure that it has power and that the power cable is plugged in. Try the Hand Control to see if it works.
- Unplug the power cable from the wall outlet. Remove dust cover. Make sure that all male and female connectors are firmly pushed into the control box. Even a slight misconnection can cause a possible malfunction. Plug the unit back into the wall and check the controls again. If one or more of the controls still do not work, contact Customer Service at 800-558-8850 to arrange replacement.



### CAUTION

Before doing any work on the unit, disconnect the power cord from the power outlet.



Be sure to read all cautions, warnings and instructions given in the manual to prevent injury to both operator and client.

## CARE &amp; CLEANING

## • Upholstery



Be sure to read all cautions, warnings and instructions given in the manual to prevent injury to both operator and client.

## IMPORTANT



## CAUTION

Before cleaning with any liquid cleaner be sure to un-plug the power cord from the power outlet. Do not operate the table in standing water.

## HOW TO CARE FOR YOUR UPHOLSTERY:

In order to maintain the quality and ensure the life of the *Relaxation Lounge™*, here are some preventative steps that can be taken to protect the upholstery on the chair and accessories.

- Clean your upholstery fabric every day. We recommend cleaning your upholstery fabric with a mild detergent to protect it against oil stains and dirt which will cause cracking and permanent stains over time. Household products you can use include a 4:1 diluted solution of 409® cleaner, Fantastik® or any non-abrasive with a non-alcohol content. Products with a high alcohol content, when used over time, will dry and eventually crack your upholstery fabric. *Note: Damage caused by these cleansers will not be covered under the warranty.*
- Avoid extremes in temperature. Very cold temperatures can make your fabric prone to cracking. Very high temperatures can cause the fabric to stretch permanently.
- Avoid exposing your products to direct sunlight for extended periods of time. This can lead to fading and cracking of the fabric.

Cushions may be replaced for a charge, but you can lengthen the life of the upholstery fabric on your lounge by following these quick and easy suggestions.

• Base  
& Controls

## GENERAL CARE OF EQUIPMENT COMPONENTS:







- Use only a mild solution of 4:1 diluted non-alcohol cleaner such as 409®, Fantastik® or some other non-abrasive cleaner to wipe down the metal and wooden surfaces. Apply with a damp (not wet) lint-free cloth. Be sure excess liquid does not drip onto, or into, any of the electrical controls or mechanisms. Be sure that the unit is unplugged before cleaning these areas. Wipe these areas carefully with a lint free cloth.
- The electronic powering mechanism of your *Relaxation Lounge™* requires no general maintenance or lubrication.
- The wooden cabinetry of the equipment can be polished with any polish approved for wooden surfaces. Spray, or pour, the polish onto a soft cloth and be careful not to get it onto other surfaces.

PRODUCT IDENTIFICATION

- Power Lead Plugs

PRODUCT IDENTIFICATION

- Plugs

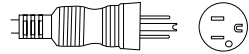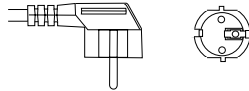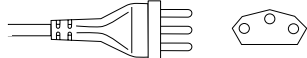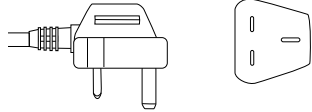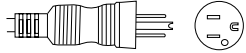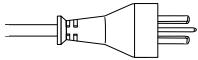Voltage (ACE)	Power Cord Plug
120v 60 Hz	3 prong grounded hospital grade; North America 
220v 50 Hz	Continental Plug; Europe/Italy 
220v 50 Hz	Swiss plug; Switzerland 
220v 50 Hz	British plug; United Kingdom 
100v 50/60 Hz	3 prong polarized plug; Japan & Korea 
220v 50 Hz	Australian Plug; Australia 

## SPECIFICATIONS



<i>Relaxation Lounge™</i>	Model No. REMPST
---------------------------	------------------

## ELECTRICAL SPECS

Designed For:	North America	Europe/Australia	Japan/Korea
Voltage/Current Service Available	120 VAC/ 20 amps	220 VAC/ 10 amps	100 VAC/ 20 amps
Input Frequency	60 Hz	50 Hz	50/60 Hz
Current Consumption w/o internal Outlets	1.8 amps	1.0 amps	1.0 amps
Power Consumption	216 watts	216 watts	216 watts
Table Actuator Voltage	24 VDC	24 VDC	24 VDC
Electric Shock Protection	Class I Equipment	Class I Equipment	Class I Equipment
Mode of Operation	Continuous Operation @ 10% Duty Cycle	Continuous Operation @ 10% Duty Cycle	Continuous Operation @ 10% Duty Cycle

## TABLE SPECS

Length:	78"	198 cm	198 cm
Width:	32"	81 cm	81 cm
Height Range*	19.5"	49 cm	49 cm
Weight*	~ 250 lbs.	~ 115 kg	~ 115 kg
Shipping Weight	~ 300 lbs.	~ 136 kg	~ 136 kg
Shipping Crate Dimensions	81" L x 26" H x 36" W	206cm x 66cm x 91cm	206cm x 66cm x 91cm
Maximum Working Load Rating	500 lbs.	227.3 kg	227.3 kg

Storage & Transport	 Temperature: 32° F - 100° F 0° C - 38° C	Humidity: 60% relative humidity During transport, DO NOT stack containers. 	Pressure: no limitations known
---------------------	---	---	--------------------------------

\* measurements are determined from floor to top of seat pad using standard padding; these measurements can be affected by natural material variations in padding.

## WARRANTY STATEMENT

<i>Relaxation Lounge</i> ™		Model No. REMPST		
FRAME/ UNDERSTRUCTURE	UPHOLSTERY /PADDING	ELECTRICAL & MOVING PARTS	RETURN SHIPPING	REPLACEMENT SHIPPING
Limited Lifetime Warranty — original owner; International : 10 years	5 years	3 year parts 1 year labor	not included	Included in Continental US; not included on International

All warranties are limited to factory-provided replacement parts, factory repair or replacement, at the discretion of OAKWORKS®.

All warranties exclude damage caused by improper set-up, accident, abuse, misuse, neglect, use for other than intended purpose or reasonable wear such as tears in the upholstery.

All warranties are invalidated by non-factory modifications and unauthorized repairs, which will immediately terminate all liability by OAKWORKS® for the product or damages caused by its use.

The buyer and its customers shall be responsible for proper set-up and use of the products as well as any supervision required for safety. In no event shall OAKWORKS® be liable for any special, indirect, consequential, incidental, exemplary or punitive damages or costs. Use of non-approved cleaning solutions voids the guarantee on all fabrics.

The warranties set forth herein are the sole and exclusive warranties provided by OAKWORKS®. There are no other warranties, representations or guarantees provided by OAKWORKS® either expressed or implied, including warranties of merchantability and fitness for a particular purpose.

Warranty valid with proof of purchase.

## LIST OF PARTS

*Relaxation Lounge™*  
MODEL REMPST

	Description
7009	Actuator LA27(2)
4046-1	Control Box (1) 110V
4046-2	Control Box (1) 220V
2092	Hand Set - 4 Button
5097	Power Cord North America
5097-02	Power Cord United Kingdom
5097-03	Power Cord Europe

Accessory  
Options

Part No.	Description
9783	Sound System
18900	Reading Light

## CONTACT INFORMATION:

### **OAKWORKS® Inc.**

923 East Wellspring Road  
New Freedom, PA 17349

Toll Free (USA only): 800-916-4613  
Phone: 717-235-6807

FAX: 717-235-6798

**[www.oakworks.com](http://www.oakworks.com)**

Manual Part # MMMNST0007 rev. 05.13.08  
Edition 1, English Sept 2007  
Edition 2, English May 2008

Printed in U.S.A.

