



## Skin Inspection Unit

### Parts/Accessories

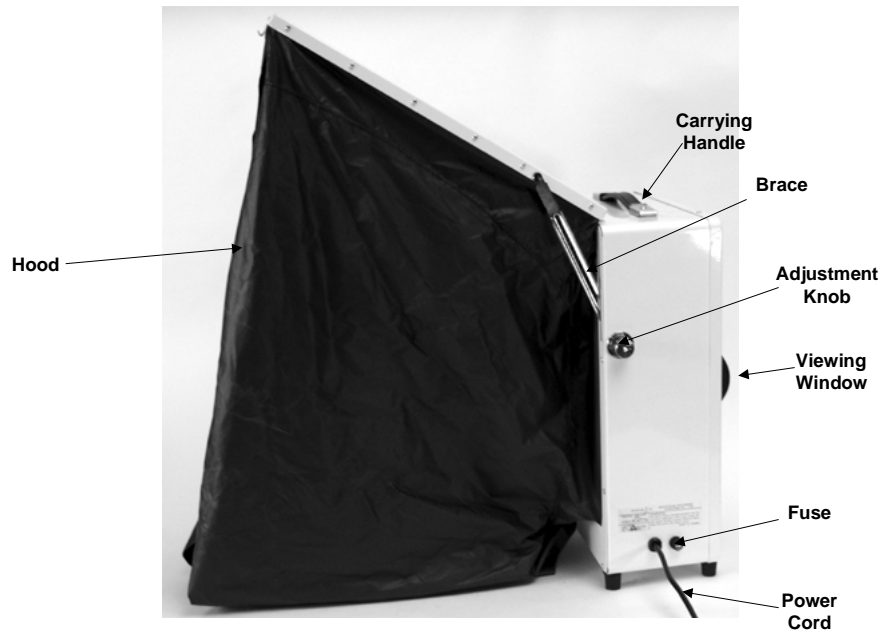

### Internal View

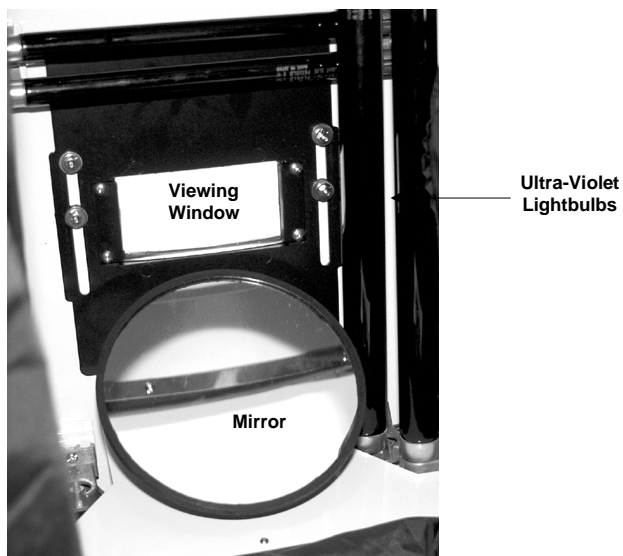

## Instructions for Use

### Setting Up Your Unit

1. Position the Skin Inspection unit on a stable, flat surface. You will want the back of the unit to be in a position where clients can sit comfortably with the curtain over their face and shoulders.
2. Release the lock located on the bottom of the Skin Inspection unit and open the lid to fully release the curtain.
3. Locate the round knob located on the side of the Skin Inspection unit. Raise the lid until it is fully extended and turn that knob counterclockwise until it locks in place.

## Specifications

### Dimensions

- 12 ¼" wide
- 17 ¼" high
- 9" deep (21 ½" deep opened up)

### Electrical

- 110 Volts
- 60 Hz
- 5 amp short fuse